



Enabling Legendary Discovery™

LEGENDplex™

Cat. No. 740492

**Mouse Immunoglobulin Isotyping Panel
(6-plex) with Filter Plate**

Cat. No. 740493

**Mouse Immunoglobulin Isotyping Panel
(6-plex) with V-bottom Plate**

Please read the entire manual before running the assay.

BioLegend.com

It is highly recommended that this manual be read in its entirety before using this product. Do not use this kit beyond the expiration date.

For Research Purposes Only. Not for use in diagnostic or therapeutic procedures. Purchase does not include or carry the right to resell or transfer this product either as a stand-alone product or as a component of another product. Any use of this product other than the permitted use without the express written authorization of BioLegend is strictly prohibited.

Table of Contents	Page
Chapter 1: KIT DESCRIPTION	3
Introduction.....	3
Principle of the Assay.....	3
Beads Usage.....	4
Storage Information.....	6
Materials Supplied.....	6
Materials to be Provided by the End-User.....	7
Precautions.....	8
Chapter 2: ASSAY PREPARATION	9
Sample Collection and Handling.....	9
Reagent Preparation.....	9
Standard Preparation.....	10
Sample Dilution.....	10
Chapter 3: ASSAY PROCEDURE	11
Performing the Assay Using a Filter Plate.....	11
Performing the Assay Using a V- bottom Plate.....	14
Chapter 4: FLOW CYTOMETER SETUP	17
Chapter 5: DATA ACQUISITION AND ANALYSIS	17
Data Acquisition.....	17
Data Analysis.....	18
Chapter 6: ASSAY CHARACTERIZATION	19
Representive Standard Curve.....	19
Assay Sensitivity.....	19
Cross-Reactivity.....	19
Linearity of Dilution.....	20
Spike Recovery.....	20

LEGENDplex™ Mouse Immunoglobulin Isotyping Panel

Intra-Assay Precision.....	21
Inter-Assay Precision.....	21
Biological Samples.....	22
TROUBLESHOOTING.....	24
PLATE MAP.....	29

Chapter 1: KIT DESCRIPTION

Introduction

Immunoglobulins, also known as antibodies, are glycoproteins produced by plasma cells. They play important roles in immune responses by recognizing, binding and neutralizing specific antigens such as bacteria, viruses and toxins. In addition, immunoglobulins can bind and mark the pathogens, a process called opsonization through which phagocyte can then engulf and eliminate the pathogens.

During immune response, plasma cells can switch from producing one immunoglobulin class to another by changing the amino acid in the constant region of the heavy chain. There are five different classes of immunoglobulins: IgM, IgG, IgE, IgA and IgD. IgG and IgA can be further divided into different subclasses based on the difference in the number of disulfide bonds and the length and flexibility of the hinge region. An increase or decrease of the immunoglobulin classes and subclasses are associated with different diseases such as lymphoproliferation or immunoglobulin-deficiency disorders. Moreover, quantitation of these classes and subclasses during and after vaccination could provide very useful information about primary and secondary immune responses to a vaccine.

The Mouse Immunoglobulin Isotyping panel is a multiplex bead-based assay panel, using fluorescence-encoded beads suitable for use on various flow cytometers. This panel allows simultaneous quantification of 6 mouse immunoglobulins, including IgG1, IgG2a, IgG2b, IgG3, IgA, and IgM. This assay panel provides broader dynamic ranges than traditional ELISA methods. The panel has been validated for use with serum, plasma and cell culture supernatant samples.

The Mouse Immunoglobulin Isotyping Panel is designed to allow flexible customization within the panel. Please visit www.biolegend.com/legendplex for more information on how to mix and match within the panel.

Principle of the Assay

BioLegend's LEGENDplex™ assays are bead-based immunoassays using the same basic principle as sandwich immunoassays.

Beads are differentiated by size and internal fluorescence intensities. Each bead set is conjugated with a specific antibody on its surface and serves as the capture beads for that particular analyte. When a selected panel of capture beads is mixed and incubated with a sample containing target analytes specific to the capture antibodies, each analyte will bind to its specific capture beads. After washing, a biotinylated detection antibody cocktail is added, and each detection antibody in the cocktail will bind to its specific analyte bound on the capture beads, thus forming capture bead-analyte-detection antibody sandwiches.

LEGENDplex™ Mouse Immunoglobulin Isotyping Panel

Streptavidin-phycoerythrin (SA-PE) is subsequently added, which will bind to the biotinylated detection antibodies, providing fluorescent signal intensities in proportion to the amount of bound analytes.


Since the beads are differentiated by size and internal fluorescence intensity on a flow cytometer, analyte-specific populations can be segregated and PE fluorescent signal quantified. The concentration of a particular analyte is determined using a standard curve generated in the same assay.

Beads Usage

The LEGENDplex bead-based assay system uses two sets of beads. Each set has a unique size that can be identified based on their forward scatter (FSC) and side scatter (SSC) profiles (Beads A and Beads B, Figure 1). Each bead set can be further resolved based on their internal fluorescence intensities. The internal dye can be detected using FL3, FL4, or APC channel, depending on the type of flow cytometer used. The smaller Beads (A) consist of 6 (shown in Figure 2) bead populations and 3 of them (A4, A5, and A7) are used for the mouse immunoglobulin isotyping panel. The larger Beads (B) consist of 7 (shown in Figure 3) bead populations and 3 of them (B2, B3 and B4) are used for this panel.

Using a total of 6 bead populations distinguished by size and internal fluorescent dye, the Mouse Immunoglobulin Isotyping Panel allows simultaneous detection of 6 immunoglobulins in a single sample. Each analyte is associated with a particular bead set as indicated (Figures 2-3 and Table 1).

Figure 1. Beads Differentiated by Size

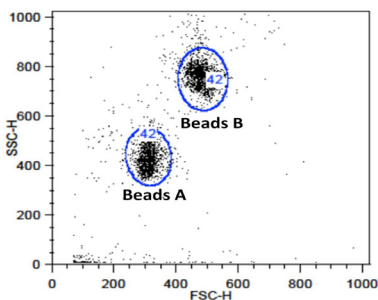

Beads A = smaller beads

Beads B = larger beads

LEGENDplex™ Mouse Immunoglobulin Isotyping Panel

Figure 2. Beads A Classification by FL4


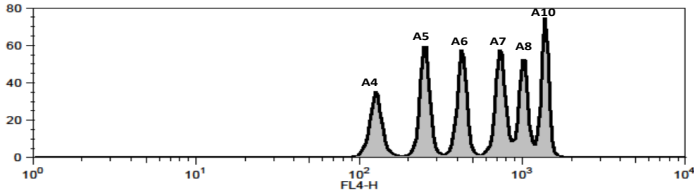


Figure 3. Beads B Classification by FL4

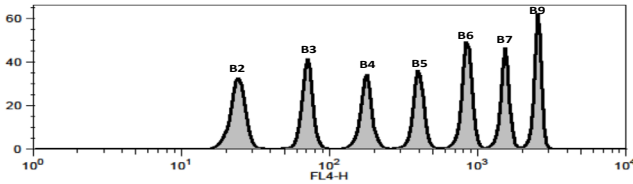

For Beads usage in the panel, please refer to Table 1 below:

Table 1. Panel Targets and Bead ID*

Target	Bead ID	Top Standard Concentrations
IgG1	A4	Note: The top standard concentrations of analytes in this panel were set at various concentrations, but may be subject to change from lot to lot (please visit biolegend.com/en-us/legendplex to download a lot-specific certificate of analysis).
IgG2a	A5	
IgG2b	B4	
IgG3	A7	
IgA	B2	
IgM	B3	

*Bead ID is used to associate a bead population to a particular analyte when using the LEGENDplex™ data analysis software program. For further information regarding the use of the program please visit biolegend.com/en-us/legendplex.

LEGENDplex™ Mouse Immunoglobulin Isotyping Panel

Storage Information

Recommended storage for all original kit components is between 2°C and 8°C. DO NOT FREEZE Pre-mixed Beads, Detection Antibodies or SA-PE.

- Once the standards have been reconstituted, immediately transfer contents into polypropylene vials. DO NOT STORE RECONSTITUTED STANDARDS IN GLASS VIALS.
- Upon reconstitution, leftover standard should be stored at $\leq -70^{\circ}\text{C}$ for use within one month. Avoid multiple (>2) freeze-thaw cycles. Discard any leftover diluted standards.

Materials Supplied

The LEGENDplex™ kit contains reagents for 100 tests, listed in the table below. When assayed in duplicate, this is enough for an 8-point standard curve and 40 samples.

Kit Components	Quantity	Volume	Part #
Setup Beads 1: FITC Beads	1 vial	1 mL	77840
Setup Beads 2: PE Beads	1 vial	1 mL	77842
Setup Beads 3: Raw Beads	1 vial	2 mL	77844
Mouse Immunoglobulin Isotyping-Panel Premixed Beads	1 bottle	3.5 mL	75065
Mouse Immunoglobulin Isotyping-Panel Detection Antibodies	1 bottle	3.5 mL	75073
Mouse Immunoglobulin Isotyping-Panel Standard Cocktail, Lyophilized	1 vial	lyophilized	75075
LEGENDplex™ SA-PE	1 bottle	3.5 mL	77743
LEGENDplex™ Assay Buffer	2 bottles	25 mL each	77562
LEGENDplex™ Wash Buffer, 20X	1 bottle	25 mL	77564
Filter Plate* or V-bottom Plate**	1 plate		76187* or 76883**
Plate Sealers	4 sheets		78101

* For assay with filter plate. ** For assay with V-bottom plate.

LEGENDplex™ Mouse Immunoglobulin Isotyping Panel

Materials to be Provided by the End-User

- A flow cytometer equipped with two lasers (e.g., a 488 nm blue laser or 532 nm green laser and a 633-635 nm red laser) capable of distinguishing 575 nm and 660 nm or a flow cytometer equipped with one laser (e.g., 488 nm blue laser) capable of distinguishing 575 nm and 670 nm.

Partial list of compatible flow cytometers:

Flow Cytometer	Reporter Channel	Channel Emission	Classification Channel	Channel Emission	Compensation needed?
BD FACSCalibur™ (single laser)	FL2	575 nm	FL3	670 nm	Yes
BD FACSCalibur™ (dual laser)	FL2	575 nm	FL4	660 nm	No*
BD Accuri™ C6	FL2	585 nm	FL4	675 nm	No*
BD FACSAria™	Yellow	575 nm	Red	660 nm	No*
BD FACSCanto™ BD FACSCanto™ II	PE	575 nm	APC	660 nm	No*
BD™ LSR, LSR II BD LSRFortessa™	PE	575-585 nm	APC	660 nm	No*
BD FACSAria™	PE	575 nm	APC	660 nm	No*
Beckman Coulter-CytoFLEX	PE	585 nm	APC	660 nm	No*

***Compensation is not required for the specified flow cytometers when set up properly.**

For setting up various flow cytometers, please visit: [www.biolegend.com/legendplex](http://www.biolegend.com/) and click on the **Instrument Setup** tab.

- Multichannel pipettes capable of dispensing 5 µL to 200 µL
- Reagent reservoirs for multichannel pipette
- Polypropylene microfuge tubes (1.5 mL)
- Laboratory vortex mixer
- Sonicator bath (e.g., Branson Ultrasonic Cleaner model #B200, or equivalent)
- Aluminum foil
- Absorbent pads or paper towels
- Plate shaker (e.g., Lab-Line Instruments model #4625, or equivalent)
- Tabletop centrifuges (e.g., Eppendorf centrifuge 5415C, or equivalent)

LEGENDplex™ Mouse Immunoglobulin Isotyping Panel

If the assay is performed in a filter plate (recommended):

- A vacuum filtration unit (Millipore MultiScreen® HTS Vacuum Manifold, cat # MSVMHTS00 or equivalent). Instructions on how to use the vacuum manifold can be found at the supplier's website.
- A vacuum source (mini vacuum pump or line vacuum, e.g., Millipore Vacuum Pump, catalog # WP6111560, or equivalent)
- If needed, additional Filter plate can be ordered from BioLegend (Cat# 740377 or 740378)

If the assay is performed in a V-bottom plate (optional):

- Centrifuge with a swinging bucket adaptor for microtiter plates (e.g., Beckman Coulter Allegra™ 6R Centrifuge with MICROPLUS CARRIER adaptor for GH3.8 and JS4.3 Rotors) .
- If needed, additional V-bottom plate can be ordered from BioLegend (Cat# 740379)

Precautions

- All blood components and biological materials should be handled as potentially hazardous. Follow universal precautions as established by the Center for Disease Control and Prevention and by the Occupational Safety and Health Administration when handling and disposing of infectious agents.
- Sodium azide has been added to some reagents as a preservative. Although the concentrations are low, sodium azide may react with lead and copper plumbing to form highly explosive metal azides. On disposal, flush with a large volume of water to prevent azide build-up.
- Do not mix or substitute reagents from different kits or lots. Reagents from different manufacturers should not be used with this kit.
- Do not use this kit beyond its expiration date.
- SA-PE and beads are light-sensitive. Minimize light exposure.

Chapter 2: ASSAY PREPARATION

Sample Collection and Handling

Preparation of Serum Samples:

- Allow the blood to clot for at least 30 minutes and centrifuge for 10 minutes at 1,000 x *g*.
- Remove serum and assay immediately or aliquot and store samples at $\leq -20^{\circ}\text{C}$. Avoid multiple (>2) freeze/thaw cycles.
- When using frozen samples, it is recommended that samples are thawed completely, mixed and centrifuged to remove particulates prior to use.

Preparation of Plasma Samples:

- Plasma collection using either EDTA, Heparin, or Citrate as an anti-coagulant are acceptable sample types to be tested with this kit. Centrifuge for 10 minutes at 1,000 x *g* within 30 minutes of blood collection.
- Remove plasma and assay immediately, or aliquot and store samples at $\leq -20^{\circ}\text{C}$. Avoid multiple (>2) freeze/thaw cycles.
- When using frozen samples, it is recommended that samples are thawed completely, mixed well and centrifuged to remove particulates.

Preparation of Cell Culture Supernatant:

- Centrifuge the sample to remove debris and assay immediately. If not possible, aliquot and store samples at $\leq -20^{\circ}\text{C}$. Avoid multiple (>2) freeze/thaw cycles.

Reagent Preparation

Preparation of Antibody-Immobilized Beads

Sonicate Pre-mixed Beads bottle for 1 minute in a sonicator bath and then vortex for 30 seconds prior to use. If no sonicator bath is available, increase the vortexing time to 1 minute to completely resuspend the beads.

Preparation of Wash Buffer

- Bring the 20X Wash Buffer to room temperature and mix to bring all salts into solution.
- Dilute 25 mL of 20X Wash Buffer with 475 mL deionized water. Store unused portions between 2°C and 8°C for up to one month.

Standard Preparation

1. Prior to use, reconstitute the lyophilized Mouse Immunoglobulin Isotyping-Panel Standard Cocktail with 250 μ L Assay Buffer.
2. Mix and allow the vial to sit at room temperature for 10 minutes, and then transfer the standard to an appropriately labeled polypropylene microfuge tube. This will be used as the top standard C7.

Note: The top standard concentrations of analytes in this panel were set at various concentrations, but may be subject to change from lot to lot (please visit biolegend.com/en-us/legendplex to download a lot-specific certificate of analysis).

3. Label 6 polypropylene microfuge tubes as C6, C5, C4, C3, C2 and C1, respectively.
4. Add 75 μ L of Assay Buffer to each of the six tubes. Prepare 1:4 dilution of the top standard by transferring 25 μ L of the top standard C7 to the C6 tube and mix well. This will be the C6 standard.
5. In the same manner, perform serial 1:4 dilutions to obtain C5, C4, C3, C2 and C1 standards (see the table below using the top standard at 6400 ng/mL as an example). Assay Buffer will be used as the 0 ng/mL standard (C0).

Tube/Standard ID	Serial Dilution	Assay Buffer to add (μ L)	Standard to add	Final Conc. (ng/mL)
C7	--	--	--	6400
C6	1:4	75	25 μ L of C7	1600
C5	1:16	75	25 μ L of C6	400
C4	1:64	75	25 μ L of C5	100
C3	1:256	75	25 μ L of C4	25
C2	1:1024	75	25 μ L of C3	6.25
C1	1:4096	75	25 μ L of C2	1.56
C0	--	75	--	0

Sample Dilution

- Serum or plasma samples should be diluted 50,000 fold with Assay Buffer as described in the table below. If further dilution is needed, the dilution should be done with Assay Buffer.

LEGENDplex™ Mouse Immunoglobulin Isotyping Panel

Sample	1 st dilution (1:200)	2 nd dilution (1:250)	Final dilution fold
Serum, plasma	2 µL + 398 µL buffer	2 µL 1 st dilution + 498 µL buffer	50,000

- For cell culture supernatant samples, the levels of analyte can vary greatly from sample to sample. While the samples can be tested without dilutions, a preliminary experiment may be required to determine the appropriate dilution factor.
- If sample dilution is desired, dilution should be done with corresponding fresh cell culture medium or Assay Buffer to ensure accurate measurement.

Chapter 3: ASSAY PROCEDURE

The LEGENDplex™ assay can be performed in a filter plate, or in a V-bottom plate.

- The in-filter plate assay procedure is recommended due to its good sample-to-sample consistency, assay robustness and ease of handling. This procedure requires a vacuum filtration unit for washing (see **Materials to be Provided by the End-User, page 7**). If you have performed bead-based multiplex assays before, your lab may already have the vacuum filtration unit set up.
- If the in-filter plate assay procedure is not possible or if you prefer, the assay can be performed in a V-bottom plate.

Performing the Assay Using a Filter Plate

- Allow all reagents to warm to room temperature (20-25°C) before use.
- Set the filter plate on an inverted plate cover at all times during assay setup and incubation steps, so that the bottom of the plate does not touch any surface. Touching a surface may cause leakage.
- Keep the plate upright during the entire assay procedure, including the washing steps, to avoid losing beads.
- The plate should be placed in the dark or wrapped with aluminum foil for all incubation steps.
- **Standards and samples should be run in duplicate and arranged on the plate in a vertical configuration convenient for data acquisition and analysis (as shown in attached PLATE MAP, page 29). Be sure to load standards in the first two columns. If an automation device is used for reading, the**

LEGENDplex™ Mouse Immunoglobulin Isotyping Panel
orientation and reading sequence should be carefully planned.

1. Pre-wet the plate by adding 100 μ L of LEGENDplex™ 1X Wash Buffer to each well and let it sit for 1 minute at room temperature. To remove the excess volume, place the plate on the vacuum manifold and apply vacuum. Do not exceed 10 inHg of vacuum. Vacuum until wells are drained (5-10 seconds). Blot excess Wash Buffer from the bottom of the plate with an absorbent pad or paper towels. Place the plate on top of the inverted plate cover.
2. Load all the wells including standards and samples wells with 25 μ L of Assay Buffer.
3. Load 25 μ L of either prepared standards to corresponding standard wells or diluted samples to sample wells.
4. Vortex mixed beads bottle for 30 seconds. Add 25 μ L of mixed beads to each well. The volume should be 75 μ L in each well after beads addition. (Note: During addition of the beads, shake mixed beads bottle intermittently to avoid bead settling).
5. Seal the plate with a plate sealer. **To avoid plate leaking, do not apply positive pressure to the sealer when sealing the plate.** Wrap the entire plate, including the inverted plate cover, with aluminum foil. Place the plate on a plate shaker, secure it and shake at approximate 500 rpm for 2 hours at room temperature.
6. **Do not invert the plate!** Place the plate on the vacuum manifold and apply vacuum as before in Step 1. Add 200 μ L of 1X Wash Buffer to each well. Remove Wash Buffer by vacuum filtration. Blot excess Wash Buffer from the bottom of the plate with an absorbent pad or paper towels. Repeat this washing step once more.
7. Add 25 μ L of Detection Antibodies to each well.
8. Seal the plate with a fresh plate sealer. Wrap the entire plate, including the inverted plate cover, with aluminum foil. Place the plate on a plate shaker and shake at approximately 500 rpm for 1 hour at room temperature.
9. **Do not vacuum!** Add 25 μ L of SA-PE to each well directly.
10. Seal the plate with a fresh plate sealer. Wrap the entire plate, including the inverted plate cover, with aluminum foil. Place the plate on a plate shaker and shake at approximate 500 rpm for 30 minutes at room temperature.
11. Repeat step 6 above.

LEGENDplex™ Mouse Immunoglobulin Isotyping Panel

12. Add 150 μL of 1X Wash Buffer to each well. Resuspend the beads on a plate shaker for 1 minute.
13. Read samples on a flow cytometer, preferably within the same day of the assay (Note: Prolonged sample storage can lead to reduced signal).


If the flow cytometer is equipped with an autosampler, read the plate directly using the autosampler. **Please be sure to program the autosampler to resuspend beads in the well immediately before taking samples. The probe height may need to be adjusted when using an autosampler.**

If an autosampler is not available, the samples can be transferred from the filter plate to micro FACS (or FACS) tubes and read manually.

Assay Procedure Summary for Filter Plate

Add 100 μL 1X Wash Buffer to filter plate wells

Vacuum to remove excess buffer

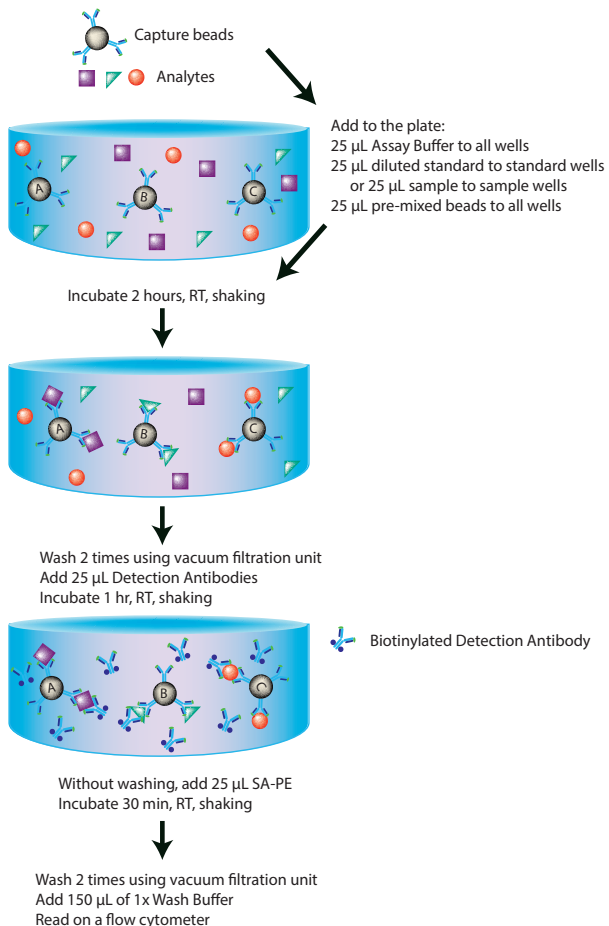

Performing the Assay Using a V-bottom Plate

- Allow all reagents to warm to room temperature (20-25°C) before use.
 - Keep the plate upright during the entire assay procedure, including the washing steps, to avoid losing beads.
 - The plate should be placed in the dark or wrapped with aluminum foil for all incubation steps.
 - Standards and samples should be run in duplicate and arranged on the plate in a vertical configuration convenient for data acquisition and analysis (as shown in attached PLATE MAP, page 29). Be sure to load standards in the first two columns of plate. If an automation device is used for reading, the orientation and reading sequence should be carefully planned.
1. Load all the wells including standards and samples wells with 25 µL of Assay Buffer.
 2. Load 25 µL of either prepared standards to corresponding standard wells or diluted samples to sample wells.
 3. Vortex mixed beads for 30 seconds. Add 25 µL of mixed beads to each well. The total volume should be 75 µL in each well after beads addition. (Note: During beads addition, shake mixed beads bottle intermittently to avoid bead settling).
 4. Seal the plate with a plate sealer. Cover the entire plate with aluminum foil to protect the plate from light. Shake at 800 rpm on a plate shaker for 2 hours at room temperature (**Depending on the shaker, the speed may need to be adjusted. The optimal speed is one that is high enough to keep beads in suspension during incubation, but not too high so it causes spill from the wells**).
 5. Centrifuge the plate at 1050 rpm (~250 g) for 5 minutes, using a swinging bucket rotor (G.H 3.8) with microplate adaptor (Please refer to **Materials to be Provided by the End-User, page 7**). Do not use excessive centrifugation speed as it may make it harder to resuspend beads in later steps. . **Make sure the timer of the centrifuge works properly and standby to make sure the centrifuge reaches preset speed.**
 6. Immediately after centrifugation, dump the supernatant into a sink by quickly inverting and flicking the plate **in one continuous and forcefull motion**. Do not worry about losing beads even if the pellet is not visible. The beads will stay in the tip of the well nicely. Blot the plate on a stack of clean paper towel and drain the remaining liquid from the well as much as possible. Be careful not to disturb the bead pellet.

Alternatively, removal of the supernatant may be completed using a multichannel pipette set at 75 µL. Try to remove as much liquid as possible

LEGENDplex™ Mouse Immunoglobulin Isotyping Panel

without removing any beads. Be sure to change pipette tips between each row or column.


7. Wash the plate by dispensing 200 μ L of washing buffer into each well and incubate for one minute. Repeat step 5 and 6 above. A second wash is optional, but may help reduce background.
8. Add 25 μ L Detection Antibodies to each well.
9. Seal the plate with a new plate sealer. Cover the entire plate with aluminum foil to protect the plate from light. Shake at 800 rpm on a plate shaker for 1 hour at room temperature.
- 10. Do not wash the plate!** Add 25 μ L of SA-PE to each well directly.
11. Seal the plate with a new plate sealer. Wrap the entire plate with aluminum foil and shake the plate on a plate shaker at approximate 800 rpm for 30 minutes at room temperature.
12. Repeat step 5 and 6
13. Wash the plate by dispensing 200 μ L of 1X Wash Buffer into each well and incubate for one minute. Repeat step 5 and 6 above. This washing step is optional but it helps to reduce the background.
14. Add 150 μ L of 1X Wash Buffer to each well. Resuspend the beads by pipetting.
15. Read samples on a flow cytometer, preferably within the same day of the assay (Note: Prolonged sample storage can lead to reduced signal).

If the flow cytometer is equipped with an autosampler, the samples can be read directly. **Please be sure to program the autosampler to resuspend beads in the well immediately before taking samples. The probe height may need to be adjusted when using an autosampler.**

If an autosampler is not available, the samples can be transferred from the plate to micro FACS (or FACS) tubes and read manually.

LEGENDplex™ Mouse Immunoglobulin Isotyping Panel

Assay Procedure Summary for V-bottom Plate

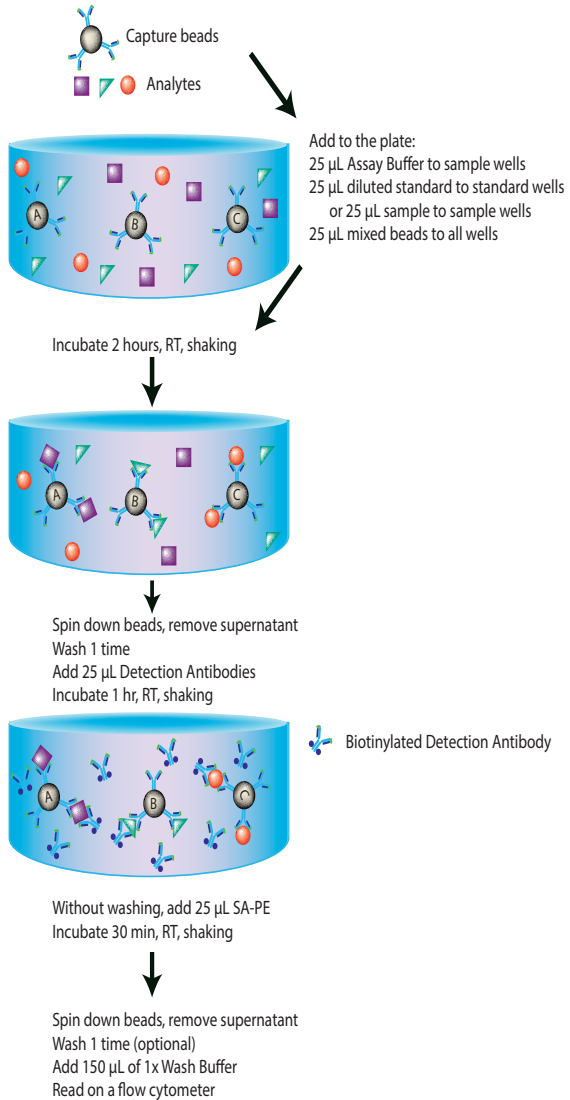

Chapter 4: FLOW CYTOMETER SETUP

In order to generate accurate data, the flow cytometer must be set up properly before data acquisition.

The setup instructions have been removed from this manual and uploaded onto our website to save paper.

To access the setup instructions, please visit: www.biolegend.com/legendplex and click on the **Instrument Setup** tab.

Chapter 5: DATA ACQUISITION AND ANALYSIS

Data Acquisition

1. Before reading samples, make sure that the flow cytometer is set up properly.
2. Create a new template or open an existing template (for details on how to create a cytometer-specific template, please refer to the Flow Cytometer Setup Guide).
3. Vortex each sample for 5 seconds before analysis.
4. Set the flow rate to low. Set the number of beads to be acquired to about 300 per analyte (e.g., acquire 1,800 beads for a 6-plex assay). Do not set to acquire total events as samples may contain large amounts of debris. Instead, create a large gate to include both Beads A and Beads B (gate A+B) and set to acquire the number of events in gate A + B. This will exclude majority of debris.

Note: Do not acquire too few or too many beads. Too few beads acquired may result in high CVs and too many beads acquired may result in slow data analysis later.

5. Read samples.

When reading samples, set the flow cytometer to setup mode first and wait until bead population is stabilized before recording or switching to acquisition mode.

To simplify data analysis using the LEGENDplex™ Data Analysis Software, read samples in the same order as shown on the PLATE MAP attached at the end of the manual. For an in-plate assay, read column by column (A1, B1, C1...A2, B2, C2...).

LEGENDplex™ Mouse Immunoglobulin Isotyping Panel

When naming data files, try to use simple names with a consecutive numbering for easy data analysis (e.g. for standards, C0.001, C0.002, C1.003, C1.004, C2.005, C2.006, C3.007, C3.008, ... C7.015, C7.016; for samples, S1.017, S1.018, S2.019, S2.020, S3.021, S3.022...)

Store all FCS files in the same folder for each assay. If running multiple assays, create a separate folder for each assay.

6. Proceed to data analysis using LEGENDplex™ Data Analysis Software when data acquisition is completed.


Data Analysis

- The assay FCS files should be analyzed using BioLegend's LEGENDplex™ data analysis software. The program is offered free of charge with the purchase of any LEGENDplex™ assay. For further information regarding access to, and use of the program please visit [biolegend.com/en-us/legendplex](https://www.biolegend.com/en-us/legendplex).

Chapter 6: ASSAY CHARACTERIZATION

Representative Standard Curve

This standard curve was generated using the LEGENDplex™ Mouse Immunoglobulin Isotyping Panel for demonstration purpose only. A standard curve must be run with each assay.

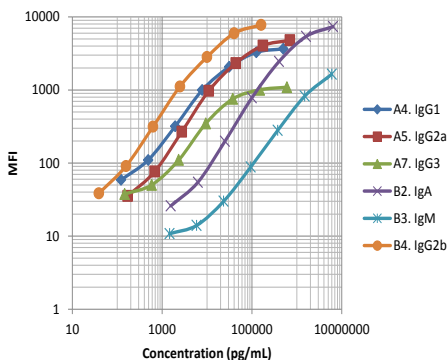

Assay Sensitivity

The assay sensitivity or minimum detectable concentration (MDC) is the theoretical limit of detection calculated using the LEGENDplex™ Data Analysis Software by applying a 5-parameter curve fitting algorithm.

Analyte	MDC (pg/mL)
IgG1	81.6 ± 64.8
IgG2a	117.0 ± 74.8
IgG2b	31.1 ± 7.2
IgG3	118.1 ± 92.0
IgA	1061.2 ± 115.6
IgM	791.2 ± 739.0

Cross-Reactivity

The following purified or recombinant mouse proteins were tested individually at 4000 ng/mL using the LEGENDplex™ Mouse Immunoglobulin

LEGENDplex™ Mouse Immunoglobulin Isotyping Panel

Isotyping Panel. There was no or negligible cross-reactivity found.

IgG1	IgG2a	IgG2b	IgG3	IgA	IgM	IgE
------	-------	-------	------	-----	-----	-----

IgG2a assay is also able to detect IgG2c. However, in C57BL/6, C57BL/10, SJL, and NOD mice, IgG2a gene is deleted and replaced with IgG2c gene. Also, IgG2b-IgG2a or IgG2b-IgG2c are the most common haplotypes in mice, and the co-existence of IgG2a and IgG2c isotypes in mice strains are rare. Therefore, the IgG2a assay will only detect either IgG2a or IgG2c in any given strain.

The kit detects all immunoglobulin classes and subclasses in mice with haplotype a and b (Igh-a and Igh-b).

Linearity of Dilution

For testing linearity of dilution, serum (n=4) and plasma samples (EDTA n=4, Heparin, Citrate, n=2 each) were first diluted 50,000 fold with Assay Buffer, respectively, then serially diluted 1:2, 1:4, 1:8 with Assay Buffer and assayed. The measured concentrations of serially diluted samples (after multiplying serial dilution factors) were then compared with the concentration of 50,000-fold diluted serum or plasma samples

Analyte	Serum	EDTA plasma	Heparin Plasma	Citrate plasma
IgG1	100%	101%	91%	86%
IgG2a	113%	108%	97%	97%
IgG2b	99%	102%	87%	77%
IgG3	102%	104%	94%	101%
IgA	98%	109%	88%	88%
IgM	103%	100%	96%	104%

Spike Recovery

For spike recovery in serum (n=4) and plasma (EDTA, Heparin and Citrate, n=4 each) samples were first diluted 50,000 fold with Assay Buffer. Then target proteins were spiked at three different levels into samples. The spiked samples were assayed, and the measured concentrations were compared with expected values.

LEGENDplex™ Mouse Immunoglobulin Isotyping Panel

Analyte	Serum	EDTA plasma	Heparin Plasma	Citrate plasma
IgG1	115%	98%	103%	90%
IgG2a	106%	95%	112%	99%
IgG2b	113%	119%	109%	118%
IgG3	81%	82%	86%	80%
IgA	92%	93%	94%	90%
IgM	91%	91%	97%	90%

Intra-Assay Precision

Two samples with different concentrations of each target protein were analyzed in one assay with 16 replicates per sample. The resulting intra-assay precision is shown below.

Analyte	Sample	Mean (ng/mL)	STDEV	%CV
IgG1	Sample 1	11.3	0.8	7
	Sample 2	2.6	0.3	12%
IgG2a	Sample 1	15.2	1.1	7%
	Sample 2	3.5	0.3	9%
IgG2b	Sample 1	3.4	0.3	7%
	Sample 2	0.8	0.1	8%
IgG3	Sample 1	13.4	1.0	7%
	Sample 2	3.3	0.4	12%
IgA	Sample 1	139.2	12.9	9%
	Sample 2	33.7	3.1	9%
IgM	Sample 1	126.5	9.1	7%
	Sample 2	29.0	3.1	11%

Inter-Assay Precision

Two samples with different concentrations of each target protein were analyzed in three independent assays with 4 replicates per sample. The resulting inter-assay precision is shown below.

LEGENDplex™ Mouse Immunoglobulin Isotyping Panel

Analyte	Sample	Mean (ng/mL)	STDEV	%CV
IgG1	Sample 1	10.5	0.8	8%
	Sample 2	2.7	0.9	7%
IgG2a	Sample 1	14.6	0.9	6%
	Sample 2	3.6	0.3	7%
IgG2b	Sample 1	3.4	0.3	8%
	Sample 2	0.9	0.04	5%
IgG3	Sample 1	12.7	0.6	5%
	Sample 2	3.4	0.3	8%
IgA	Sample 1	135.7	8.9	7%
	Sample 2	34.0	1.9	6%
IgM	Sample 1	124.4	8.6	7%
	Sample 2	30.7	2.1.	7%

Biological Samples

Serum

Normal mouse serum from four different mouse strains (C57BL/6, BALB/c, Swiss Webster and CD-1, pooled n ≥ 10 each) were tested for endogenous levels of different mouse immunoglobulin. The concentrations (mg/mL) measured are shown below.

Analyte	CD-1	Swiss Webster	BALB/c	C57BL/6
IgG1	0.8	0.3	1.0	1.1
IgG2a	1.5	0.2	1.2	2.3
IgG2b	1.0	0.07	0.3	2.1
IgG3	0.08	0.04	0.6	0.1
IgA	0.2	0.1	0.2	0.8
IgM	0.4	0.3	1.2	0.5

LEGENDplex™ Mouse Immunoglobulin Isotyping Panel

Plasma

Normal mouse EDTA plasma from four different mouse strains (C57BL/6, BALB/c, Swiss Webster and CD-1, pooled $n \geq 10$ each) were tested for endogenous levels of different mouse immunoglobulin. The concentrations (mg/mL) measured are shown below.

Analyte	CD-1	Swiss Webster	BALB/c	C57BL/6
IgG1	0.6	1.6	1.0	0.3
IgG2a	0.8	2.6	1.1	0.6
IgG2b	0.4	1.8	0.3	0.4
IgG3	0.04	0.1	0.1	0.05
IgA	0.2	0.5	0.5	0.1
IgM	0.2	0.6	0.9	0.2

Citrate Plasma

Normal mouse Citrate plasma from four different mouse strains (C57BL/6, BALB/c, Swiss Webster and CD-1, pooled $n \geq 10$ each) were tested for endogenous levels of different mouse immunoglobulin. The concentrations (mg/mL) measured are shown below.

Analyte	CD-1	Swiss Webster	BALB/c	C57BL/6
IgG1	0.8	1.2	0.4	1.6
IgG2a	1.2	3.0	0.2	3.5
IgG2b	0.8	2.4	0.08	2.3
IgG3	0.06	0.2	0.05	0.3
IgA	0.1	0.6	0.1	0.3
IgM	0.2	0.8	0.3	1.7

LEGENDplex™ Mouse Immunoglobulin Isotyping Panel

Heparin Plasma

Normal mouse Heparin plasma from four different mouse strains (C57BL/6, BALB/c, Swiss Webster and CD-1, pooled $n \geq 10$ each) were tested for endogenous levels of different mouse immunoglobulin. The concentrations (mg/mL) measured are shown below.

Analyte	CD-1	Swiss Webster	BALB/c	C57BL/6
IgG1	1.7	0.5	1.1	1.3
IgG2a	3.2	0.2	2.6	1.6
IgG2b	2.3	0.1	1.7	0.5
IgG3	0.2	0.05	0.1	0.2
IgA	0.7	0.06	0.5	0.6
IgM	0.7	0.6	0.4	1.6

TROUBLESHOOTING

Problem	Possible Cause	Solution
Bead population shifting upward or downward during acquisition	The strong PE signal from high concentration samples or standards may spill over to classification Channel (e.g., FL3/FL4/APC) and mess up the bead separation.	Optimize instrument settings using Kit Setup Beads, and make appropriate compensation between channels.
Filter plate will not vacuum or some wells clogged	Vacuum pressure is insufficient or vacuum manifold does not seal properly.	Increase vacuum pressure such that 0.2 mL buffer can be suctioned in 3-5 seconds. Clean the vacuum manifold and make sure no debris on the manifold. Press down the plate on the manifold to make a good seal.
	Samples have insoluble particles or sample is too viscous (e.g., serum and plasma samples)	<p>Centrifuge samples just prior to assay setup and use supernatant. If high lipid content is present, remove lipid layer after centrifugation. Sample may need dilution if too viscous.</p> <p>If some wells are still clogged during washing, try the following:</p> <ol style="list-style-type: none"> 1). Add buffer to all the wells, pipette up and down the clogged wells and vacuum again. 2). Use a piece of clean wipe, wipe the under side of the clogged wells and vacuum again. 3). Take a thin needle (e.g., insulin needle), while holding the plate upward, poke the little hole under each of the clogged wells and vacuum again. Do not poke too hard or too deep as it may damage the filter and cause leaking.
	Filter plate was used without pre-wet.	Pre-wet plate with wash buffer before running the assay.

LEGENDplex™ Mouse Immunoglobulin Isotyping Panel

Insufficient bead count or slow reading	Beads inappropriately prepared	Sonicate bead vials and vortex just prior to addition. Agitate mixed beads intermittently in reservoir while pipetting this into the plate.
	Samples cause beads aggregation due to particulate matter or viscosity.	Centrifuge samples just prior to assay setup and use supernatant. If high lipid content is present, remove lipid layer after centrifugation. Sample may need dilution if too viscous.
	Beads were lost during washing for in-tube assay	Make sure beads are spun down by visually check the pellet (beads are in light blue or blue color). Be very careful when removing supernatant during washing.
	Probe might be partially clogged.	Sample probe may need to be cleaned, or if needed, probe should be removed and sonicated.
Plate leaked	Vacuum pressure set too high	Adjust vacuum pressure such that 0.2 mL buffer can be suctioned in 3-5 seconds. Do not exceed 10" Hg of vacuum.
	Plate set directly on table or absorbent towels during incubations or reagent additions	Set plate on plate holder or raised edge so bottom of filter is not touching any surface.
	Liquid present on the under side of the plate after vacuum	After washing, press down plate firmly on a stack of clean paper towels to dry the underside of the plate.
	Pipette touching and damaged plate filter during additions.	Pipette to the side of wells.
High Background	Background wells were contaminated	Avoid cross-well contamination by changing tips between pipetting when performing the assay using a multichannel pipette.
	Insufficient washes	The background may be due to non-specific binding of SA-PE. Increase number of washes.
Debris (FSC/SSC) during sample acquisition	Debris or platelet may exist in sample solution.	Centrifuge samples before analyzing samples. Remove platelet as much as possible.

LEGENDplex™ Mouse Immunoglobulin Isotyping Panel

Variation between duplicate samples	Beads aggregation	Sonicate and vortex the Beads prior to use.
	Multichannel pipette may not be calibrated or inconsistent pipetting	Calibrate Pipette. Ensure good pipetting practice. Prime pipette before use may help.
	Plate washing was not uniform	Make sure all reagents are vacuumed out completely in all wash steps.
	Samples may contain particulate matters.	Centrifuge samples just prior to assay setup and use supernatant. If high lipid content is present, remove lipid layer after centrifugation. Sample may need dilution if too viscous.
Low or poor standard curve signal	The standard was incorrectly reconstituted, stored or diluted	Follow the protocol to reconstitute, store and dilute standard. Double check your calculation.
	Wrong or short incubation time	Ensure the time of all incubations was appropriate.
Signals too high, standard curves saturated	PMT value for FL2/PE set too high	Make sure the PMT setting for the reporter channel is appropriate
	Plate incubation time was too long	Use shorter incubation time.
Sample readings are out of range	Samples contain no or below detectable levels of analyte	Make sure the experiment to generate the samples worked. Use proper positive controls.
	Samples concentrations higher than highest standard point.	Dilute samples and analyze again.
	Standard curve was saturated at higher end of curve.	Make sure the PMT setting for the reporter channel is appropriate. Use shorter incubation time if incubation time was too long
Missed beads populations during reading, or distribution is unequal	Sample may cause some beads to aggregate.	Centrifuge samples just prior to assay setup and use supernatant. If high lipid content is present, remove lipid layer after centrifugation. Sample may need dilution if too viscous.
	Beads populations are not mixed properly	Make sure all bead populations are mixed. and in similar numbers.

LEGENDplex™ Mouse Immunoglobulin Isotyping Panel

PLATE MAP (for in-plate assay)

	1	2	3	4	5	6	7	8	9	10	11	12
A	C0	C4	Sample1	Sample5	Sample 9	Sample 13	Sample 17	Sample 21	Sample 25	Sample 29	Sample 33	Sample 37
B	C0	C4	Sample1	Sample5	Sample 9	Sample 13	Sample 17	Sample 21	Sample 25	Sample 29	Sample 33	Sample 37
C	C1	C5	Sample2	Sample6	Sample 10	Sample 14	Sample 18	Sample 22	Sample 26	Sample 30	Sample 34	Sample 38
D	C1	C5	Sample2	Sample6	Sample 10	Sample 14	Sample 18	Sample 22	Sample 26	Sample 30	Sample 34	Sample 38
E	C2	C6	Sample3	Sample7	Sample 11	Sample 15	Sample 19	Sample 23	Sample 27	Sample 31	Sample 35	Sample 39
F	C2	C6	Sample3	Sample7	Sample 11	Sample 15	Sample 19	Sample 23	Sample 27	Sample 31	Sample 35	Sample 39
G	C3	C7	Sample4	Sample8	Sample 12	Sample 16	Sample 20	Sample 24	Sample 28	Sample 32	Sample 36	Sample 40
H	C3	C7	Sample4	Sample8	Sample 12	Sample 16	Sample 20	Sample 24	Sample 28	Sample 32	Sample 36	Sample 40



LEGENDplex™ Kits are manufactured by **BioLegend**

8999 BioLegend Way

San Diego, CA 92121

Tel: 1.858.768.5800

Tel: US & Canada Toll-Free: 1.877.Bio-Legend (1.877.246.5343)

Fax: 1.877.455.9587

Email: info@biolegend.com

biolegend.com

For a complete list of world-wide BioLegend offices and distributors,
please visit our website at: biolegend.com