

GMP FITC anti-human CD3 Antibody

Catalog# / Size	260028 / 100 tests
Clone	SK7
Workshop	HCDM listed
Other Names	T3, CD3ε
Isotype	Mouse IgG1, κ
Description	CD3ε is a 20 kD chain of the CD3/T-cell receptor (TCR) complex, which is composed of two CD3ε, one CD3γ, one CD3δ, one CD3ζ (CD247), and a T-cell receptor (α/β or γ/δ) heterodimer. It is found on all mature T cells, NK T cells, and some thymocytes. CD3, also known as T3, is a member of the immunoglobulin superfamily that plays a role in antigen recognition, signal transduction, and T cell activation.

Product Details

Reactivity	Human
Antibody Type	Monoclonal
Host Species	Mouse
Formulation	Phosphate-buffered solution, pH 7.2, containing 0.09% sodium azide, 0.2% (w/v) BSA (origin USA) and a stabilizer.
Preparation	The antibody was purified by affinity chromatography and conjugated with FITC under optimal conditions.
Concentration	200 µg/mL
Storage & Handling	The antibody solution should be stored undiluted between 2°C and 8°C, and protected from prolonged exposure to light. Do not freeze.
Application	FC - Quality tested
Recommended Usage	Each lot of this antibody is quality control tested by immunofluorescent staining with flow cytometric analysis . For flow cytometric staining, the suggested use of this reagent is 5 µL per million cells in 100 µL staining volume or 5 µL per 100 µL of whole blood. It is recommended that the reagent be titrated for optimal performance for each application.
Excitation Laser	Blue Laser (488 nm)
Application Notes	Additional reported application (for the relevant formats) include: immunohistochemical staining of frozen tissue sections ^{4,5,8} , immunofluorescent staining ⁶ , and Western blotting ³ .

Application References

(PubMed link indicates BioLegend citation)

1. Kan EA, et al. 1983. J. Immunol. 131:536.
2. Wood GS, et al. 1985. Am. J. Pathol. 120:371.
3. Van Dongen JJM, et al. 1988. Blood 71:603. (WB)
4. Haringman JJ, et al. 2005. Arthritis Res. Ther. 7:R862. (IHC)
5. Carbone A, et al. 1999. Blood 93:2319. (IHC)
6. Goyal JJ, et al. 2006. J. Histochem. Cytochem. 54:75. (IF)
7. Rutjens E, et al. 2007. J. Immunol. 178:1702.
8. Kap Y, et al. 2009. J. Histochem. Cytochem. 57:1159. (IHC)
9. Yoshino N, et al. 2000. Exp. Anim. (Tokyo) 49:97. (FC)

Disclaimer

GMP RUO Flow Cytometry Antibodies. BioLegend GMP RUO fluorophore conjugated antibodies are manufactured in a dedicated GMP facility and compliant with ISO 13485:2016. For research use only. Not for use in diagnostic or therapeutic procedures. Our processes include:

- Batch-to-batch consistency

- Material traceability
- Documented procedures
- Documented employee training
- Equipment maintenance and monitoring records
- Lot-specific certificates of analysis
- Quality audits per ISO 13485:2016
- QA review of released products

Antigen Details

Structure	Ig superfamily, with the subunits of CD3 γ , CD3 δ , CD3 ζ , (CD247) and TCR (α/β or γ/δ) forms CD3/TCR complex, 20 kD
Distribution	Mature T and NK T cells, during thymocyte differentiation
Function	Antigen recognition, signal transduction, T cell activation
Ligand/Receptor	Peptide antigen bound to MHC
Cell Type	NKT cells, T cells, Tregs
Biology Area	Immunology, Innate Immunity
Molecular Family	CD Molecules, TCRs
Antigen References	1. Barclay N, et al. 1993. The Leucocyte FactsBook. Academic Press. San Diego. 2. Beverly P, et al. 1981. Eur. J. Immunol. 11:329. 3. Lanier L, et al. 1986. J. Immunol. 137:2501.
Gene ID	916


Related Protocols

[Cell Surface Flow Cytometry Staining Protocol](#)

Other Formats

APC/Fire™ 750 anti-human CD3, Biotin anti-human CD3, Purified anti-human CD3, FITC anti-human CD3, PE anti-human CD3, Alexa Fluor® 488 anti-human CD3, APC anti-human CD3, PerCP/Cyanine5.5 anti-human CD3, PerCP anti-human CD3, PE/Cyanine7 anti-human CD3, APC/Cyanine7 anti-human CD3, Alexa Fluor® 700 anti-human CD3, Pacific Blue™ anti-human CD3, Alexa Fluor® 647 anti-human CD3, Brilliant Violet 510™ anti-human CD3, Brilliant Violet 421™ anti-human CD3, Brilliant Violet 605™ anti-human CD3, Brilliant Violet 711™ anti-human CD3, Brilliant Violet 785™ anti-human CD3, PE/Dazzle™ 594 anti-human CD3, Brilliant Violet 750™ anti-human CD3, TotalSeq™-A0049 anti-human CD3, TotalSeq™-C0049 anti-human CD3, Spark Blue™ 550 anti-human CD3, TotalSeq™-B0049 anti-human CD3, Alexa Fluor® 660 anti-human CD3, APC/Fire™ 810 anti-human CD3, Spark NIR™ 685 anti-human CD3, PE/Fire™ 640 anti-human CD3, PE/Fire™ 700 anti-human CD3

Product Data

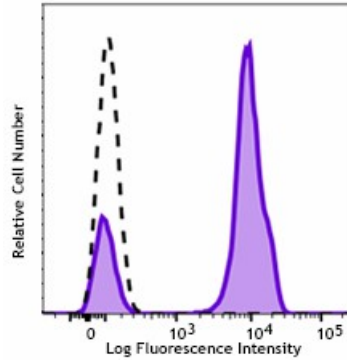

Typical results from human peripheral blood lymphocytes stained either with SK7 FITC used at 5 μ L/test (filled histogram) or with isotype control (open histogram).

Symbols Glossary*

Symbol	Meaning	Symbol Title	Symbol No.	Symbol	Meaning	Symbol Title	Symbol No.
	Catalog number	Catalogue number	5.1.6		Indicates the need for the user to consult the instructions for use.	Consult instructions for use	5.4.3
	Indicates the temperature limits to which the medical device can be safely exposed.	Temperature limit	5.3.7		Indicates a medical device that needs protection from light sources.	Keep away from sunlight	5.3.2
	Indicates the upper limit of temperature to which the medical device can be safely exposed.	Upper limit of temperature	5.3.6		Indicates the date after which the medical device is not to be used.	Use-by date	5.1.4
	Indicates the medical device manufacturer.	Manufacturer	5.1.1		Indicates the authorized representative in the European Community.	Authorized representative in the European Community	5.1.2
	Indicates the manufacturer's batch code so that the batch or lot can be identified.	Batch code	5.1.5		Indicates a medical device that is intended to be used as an in vitro diagnostic medical device.	<i>In vitro</i> diagnostic medical device	5.5.1

* Symbol information is from EN ISO 15223-1:2016 Medical devices – Symbols to be used with medical device labels, labelling and information to be supplied – Part 1: General requirements

For research use only. Not for resale. BioLegend will not be held responsible for patent infringement or other violations that may occur with the use of our products.

*These products may be covered by one or more Limited Use Label Licenses (see the BioLegend Catalog or our website, www.biolegend.com/ordering#license). BioLegend products may not be transferred to third parties, resold, modified for resale, or used to manufacture commercial products, reverse engineer functionally similar materials, or to provide a service to third parties without written approval of BioLegend. By use of these products you accept the terms and conditions of all applicable Limited Use Label Licenses.

BioLegend Inc., 8999 BioLegend Way, San Diego, CA 92121 www.biolegend.com
Toll-Free Phone: 1-877-Bio-Legend (246-5343) Phone: (858) 768-5800 Fax: (877) 455-9587