

GMP PE anti-human CD16 Antibody

Catalog# / Size	260096 / 100 tests
Clone	3G8
Workshop	V NK80
Other Names	FcγRIII, Fc gamma receptor, Fc gamma receptor 3
Isotype	Mouse IgG1, κ
Description	CD16 is known as low affinity IgG receptor III (FcγRIII). It is expressed as two distinct forms (CD16a and CD16b). CD16a (FcγRIIIA) is a 50-65 kD polypeptide-anchored transmembrane protein. It is expressed on the surface of NK cells, activated monocytes, macrophages, and placental trophoblasts in humans. CD16b (FcγRIIIB) is a 48 kD glycosylphosphatidylinositol (GPI)-anchored protein. Its extracellular domain is over 95% homologous to that of CD16a, and it is expressed specifically on neutrophils. CD16 binds aggregated IgG or IgG-antigen complex which functions in NK cell activation, phagocytosis, and antibody-dependent cell-mediated cytotoxicity (ADCC).

Product Details

Reactivity	Human
Antibody Type	Monoclonal
Host Species	Mouse
Immunogen	Human PMN cells
Formulation	Phosphate-buffered solution, pH 7.2, containing 0.09% sodium azide, 0.2% (w/v) BSA (origin USA) and a stabilizer.
Preparation	The antibody was purified by affinity chromatography and conjugated with PE under optimal conditions.
Concentration	100 µg/mL
Storage & Handling	The antibody solution should be stored undiluted between 2°C and 8°C, and protected from prolonged exposure to light. Do not freeze.
Application	FC - Quality tested
Recommended Usage	Each lot of this antibody is quality control tested by immunofluorescent staining with flow cytometric analysis . For flow cytometric staining, the suggested use of this reagent is 5 µL per million cells in 100 µL staining volume or 5 µL per 100 µL of whole blood. It is recommended that the reagent be titrated for optimal performance for each application.
Excitation Laser	Blue Laser (488 nm) Green Laser (532 nm)/Yellow-Green Laser (561 nm)
Application Notes	<p>The 3G8 antibody clone blocks neutrophil phagocytosis and stimulates NK cell proliferation. It has been reported that this clone interacts with the FcγRIIIa and FcγRIIIb receptors causing neutrophil activation and aggregation¹⁸. Due to this phenomenon staining in whole blood may cause a reduction in the number of granulocytes or alter their scatter profile.</p> <p>Additional reported applications (for the relevant formats) include: immunohistochemical staining of acetone-fixed frozen tissue sections⁶, immunoprecipitation³, stimulation of NK cell proliferation⁴, blocking of phagocytosis⁵, and blocking of immunoglobulin binding to FcγRIII^{7,8}. The Ultra-LEAF™ purified antibody (Endotoxin < 0.01 EU/µg, Azide-Free, 0.2 µm filtered) is recommended for functional assays (Cat. No. 302049, 302050, 302057, 302058).</p>
Application References	<ol style="list-style-type: none"> 1. Knapp W, <i>et al.</i> Eds. 1989. Leucocyte Typing IV. Oxford University Press. New York. 2. Schlossman S, <i>et al.</i> Eds. 1995. Leucocyte Typing V. Oxford University Press. New York. 3. Edberg J, <i>et al.</i> 1997. <i>J. Immunol.</i> 159:3849. (IP) 4. Hoshino S, <i>et al.</i> 1991. <i>Blood</i> 78:3232. (Stim)
(PubMed link indicates BioLegend citation)	

5. Tamm A, *et al.* 1996. *Immunol.* 157:1576. (Block)
6. Da Silva DM, *et al.* 2001. *Int. Immunol.* 13:633. (IHC)
7. Holl V, *et al.* 2004. *J. Immunol.* 173:6274. (Block)
8. Hober D, *et al.* 2002. *J. Gen. Virol.* 83:2169. (Block)
9. Brainard DM, *et al.* 2009. *J. Virol.* 83:7305. [PubMed](#)
10. Smed-Sørensen A, *et al.* 2008. *Blood* 111:5037. (Block) [PubMed](#)
11. Timmerman KL, *et al.* 2008. *J. Leukoc. Biol.* 84:1271. (FC) [PubMed](#)
12. Yoshino N, *et al.* 2000. *Exp. Anim. (Tokyo)* 49:97. (FC)
13. Rout N, *et al.* 2010. *PLoS One* 5:e9787. (FC)
14. Kim WK, *et al.* 2006. *Am. J. Pathol.* 168:822. (FC)
15. Boltz A, *et al.* 2011. *J. Biol Chem.* 286:21896. [PubMed](#)
16. Wu Z, *et al.* 2013. *J. Virol.* 87:7717. [PubMed](#)
17. Peterson VM, *et al.* 2017. *Nat. Biotechnol.* 35:936. (PG)
18. Vossebeld PJ, *et al.* 1997. *Biochem J.* 323:87-94 (Stim)

Disclaimer

GMP RUO Flow Cytometry Antibodies. BioLegend GMP RUO fluorophore conjugated antibodies are manufactured in a dedicated GMP facility and compliant with ISO 13485:2016. For research use only. Not for use in diagnostic or therapeutic procedures. Our processes include:

- Batch-to-batch consistency
- Material traceability
- Documented procedures
- Documented employee training
- Equipment maintenance and monitoring records
- Lot-specific certificates of analysis
- Quality audits per ISO 13485:2016
- QA review of released products

Antigen Details


Structure	Ig superfamily, transmembrane form (50-65 kD) or GPI-linked form (48 kD)
Distribution	NK cells, activated monocytes, macrophages, neutrophils
Function	Low affinity IgG Fc receptor, phagocytosis, ADCC
Ligand/Receptor	Aggregated IgG, IgG-antigen complex
Cell Type	Dendritic cells, Macrophages, Monocytes, Neutrophils, NK cells
Biology Area	Immunology, Innate Immunity
Molecular Family	CD Molecules, Fc Receptors
Antigen References	<ol style="list-style-type: none"> 1. Fleit H, <i>et al.</i> 1982. <i>P. Natl. Acad. Sci. USA</i> 79:3275. 2. Stroncek D, <i>et al.</i> 1991. <i>Blood</i> 77:1572. 3. Wirthmueller U, <i>et al.</i> 1992. <i>J. Exp. Med.</i> 175:1381.
Gene ID	2214

Related Protocols

[Cell Surface Flow Cytometry Staining Protocol](#)

Other Formats

APC anti-human CD16, Biotin anti-human CD16, FITC anti-human CD16, Brilliant Violet 711™ anti-human CD16, PE anti-human CD16, PE/Cyanine5 anti-human CD16, Purified anti-human CD16, APC/Cyanine7 anti-human CD16, PE/Cyanine7 anti-human CD16, Alexa Fluor® 488 anti-human CD16, Alexa Fluor® 647 anti-human CD16, Pacific Blue™ anti-human CD16, Alexa Fluor® 700 anti-human CD16, PerCP/Cyanine5.5 anti-human CD16, PerCP anti-human CD16, Brilliant Violet 421™ anti-human CD16, Brilliant Violet 570™ anti-human CD16, Brilliant Violet 605™ anti-human CD16, Brilliant Violet 650™ anti-human CD16, Brilliant Violet 785™ anti-human CD16, Brilliant Violet 510™ anti-human CD16, Ultra-LEAF™ Purified anti-human CD16, Purified anti-human CD16 (Maxpar® Ready), PE/Dazzle™ 594 anti-human CD16, APC/Fire™ 750 anti-human CD16, TotalSeq™-A0083 anti-human CD16, TotalSeq™-B0083 anti-human CD16, TotalSeq™-C0083 anti-human CD16, PE/Fire™ 640 anti-human CD16, Spark YG™ 581 anti-human CD16, TotalSeq™-D0083 anti-human CD16, APC/Fire™ 810 anti-human CD16, GMP APC anti-human CD16, GMP PE/Dazzle™ 594 anti-human CD16

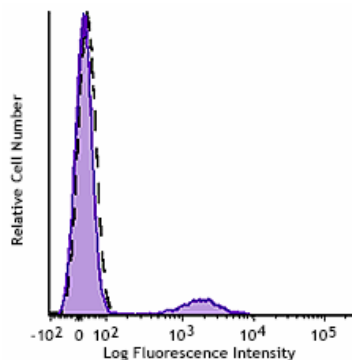

Typical results from human peripheral blood lymphocytes stained either with 3G8 PE used at 5 μ L/test (filled histogram) or with an isotype control (open histogram).

Symbols Glossary*

Symbol	Meaning	Symbol Title	Symbol No.	Symbol	Meaning	Symbol Title	Symbol No.
	Catalog number	Catalogue number	5.1.6		Indicates the need for the user to consult the instructions for use.	Consult instructions for use	5.4.3
	Indicates the temperature limits to which the medical device can be safely exposed.	Temperature limit	5.3.7		Indicates a medical device that needs protection from light sources.	Keep away from sunlight	5.3.2
	Indicates the upper limit of temperature to which the medical device can be safely exposed.	Upper limit of temperature	5.3.6		Indicates the date after which the medical device is not to be used.	Use-by date	5.1.4
	Indicates the medical device manufacturer.	Manufacturer	5.1.1		Indicates the authorized representative in the European Community.	Authorized representative in the European Community	5.1.2
	Indicates the manufacturer's batch code so that the batch or lot can be identified.	Batch code	5.1.5		Indicates a medical device that is intended to be used as an in vitro diagnostic medical device.	<i>In vitro</i> diagnostic medical device	5.5.1

* Symbol information is from EN ISO 15223-1:2016 Medical devices – Symbols to be used with medical device labels, labelling and information to be supplied – Part 1: General requirements

For research use only. Not for resale. BioLegend will not be held responsible for patent infringement or other violations that may occur with the use of our products.

*These products may be covered by one or more Limited Use Label Licenses (see the BioLegend Catalog or our website, www.biolegend.com/ordering#license). BioLegend products may not be transferred to third parties, resold, modified for resale, or used to manufacture commercial products, reverse engineer functionally similar materials, or to provide a service to third parties without written approval of BioLegend. By use of these products you accept the terms and conditions of all applicable Limited Use Label Licenses.

BioLegend Inc., 8999 BioLegend Way, San Diego, CA 92121 www.biolegend.com
 Toll-Free Phone: 1-877-Bio-Legend (246-5343) Phone: (858) 768-5800 Fax: (877) 455-9587