

GMP PE/Cyanine7 anti-human CD8 Antibody

Catalog# / Size	260106 / 100 tests
Clone	SK1
Other Names	T8, Leu2
Isotype	Mouse IgG1, κ
Description	CD8a is a 32-34 kD type I glycoprotein. It forms a homodimer (CD8a/a) or heterodimer (CD8a/b) with CD8b. CD8, also known as T8 and Leu2, is a member of the immunoglobulin superfamily found on the majority of thymocytes, a subset of peripheral blood T cells, and NK cells (which express almost exclusively CD8a homodimers). CD8 acts as a co-receptor with MHC class I-restricted T cell receptors in antigen recognition and T cell activation and has been shown to play a role in thymic differentiation. Two domains in CD8a are important for function: the extracellular IgSF domain binds the α_3 domain of MHC class I and the cytoplasmic CXCP motif binds the tyrosine kinase p56 Lck.

Product Details

Reactivity	Human
Antibody Type	Monoclonal
Host Species	Mouse
Formulation	Phosphate-buffered solution, pH 7.2, containing True-Stain Monocyte Blocker™, 0.09% sodium azide and 0.2% (w/v) BSA (origin USA).
Preparation	The antibody was purified by affinity chromatography and conjugated with PE/Cyanine7 under optimal conditions.
Concentration	100 μ g/mL
Storage & Handling	The antibody solution should be stored undiluted between 2°C and 8°C, and protected from prolonged exposure to light. Do not freeze.
Application	FC - Quality tested
Recommended Usage	Each lot of this antibody is quality control tested by immunofluorescent staining with flow cytometric analysis . For flow cytometric staining, the suggested use of this reagent is 5 μ L per million cells in 100 μ L staining volume or 5 μ L per 100 μ L of whole blood. It is recommended that the reagent be titrated for optimal performance for each application.
Excitation Laser	Blue Laser (488 nm) Green Laser (532 nm)/Yellow-Green Laser (561 nm)
Application Notes	Clone SK1 recognizes the a chain of CD8. Additional reported applications (for the relevant formats) include: proteogenomics ⁸ , immunohistochemistry of acetone-fixed frozen tissue sections, and spatial biology (IBEX) ^{9,10} . This clone was tested in-house and does not demonstrate utility for formalin-fixed paraffin-embedded (FFPE) human tonsil sections.
Application References	<ol style="list-style-type: none"> 1. Ledbetter JA, <i>et al.</i> 1981. <i>J. Exp. Med.</i> 153:310. 2. Campanelli R, <i>et al.</i> 2002. <i>Intl. Immunol.</i> 14:39. 3. Evans RL, <i>et al.</i> 1981. <i>Immunol.</i> 78:544. 4. Wooldridge L, <i>et al.</i> 2005. <i>J. Bio. Chem.</i> 280:27491. 5. Ch'el IL, <i>et al.</i> 2011. <i>J Exp Med.</i> 208:633. PubMed 6. Carbone A, <i>et al.</i> 1999. <i>Blood</i> 93:2319. (IHC-F) 7. Ahmed A, <i>et al.</i> 2001. <i>J. Pathol.</i> 193:383. (IHC) 8. Peterson VM, <i>et al.</i> 2017. <i>Nat. Biotechnol.</i> 35:936. (PG) 9. Radtke AJ, <i>et al.</i> 2020. <i>Proc Natl Acad Sci USA.</i> 117:33455-33465. (SB) PubMed 10. Radtke AJ, <i>et al.</i> 2022. <i>Nat Protoc.</i> 17:378-401. (SB) PubMed
(PubMed link indicates BioLegend citation)	
Disclaimer	GMP RUO Flow Cytometry Antibodies. BioLegend GMP RUO fluorophore conjugated antibodies are manufactured in a dedicated GMP facility and compliant with ISO 13485:2016. For research use only. Not for use in diagnostic or therapeutic procedures. Our processes include:

- Batch-to-batch consistency
- Material traceability
- Documented procedures
- Documented employee training
- Equipment maintenance and monitoring records
- Lot-specific certificates of analysis
- Quality audits per ISO 13485:2016
- QA review of released products

Antigen Details


Structure	Ig superfamily, homodimer or heterodimer with CD8b, 32-34 kD
Distribution	Majority of thymocytes, T cell subset, NK cells
Function	MHC class I co-receptor, thymic differentiation, T cell activation
Ligand/Receptor	MHC Class I molecules
Cell Type	NK cells, T cells, Thymocytes
Biology Area	Immunology
Molecular Family	CD Molecules
Antigen References	1. Barclay N, <i>et al.</i> 1993. <i>The Leucocyte Antigen FactsBook</i> . Academic Press Inc. San Diego.
Gene ID	925

Related Protocols

[Cell Surface Flow Cytometry Staining Protocol](#)

Other Formats

Alexa Fluor® 647 anti-human CD8, Brilliant Violet 650™ anti-human CD8, Purified anti-human CD8, FITC anti-human CD8, PE anti-human CD8, PerCP anti-human CD8, PerCP/Cyanine5.5 anti-human CD8, PE/Cyanine7 anti-human CD8, APC/Cyanine7 anti-human CD8, Alexa Fluor® 488 anti-human CD8, Pacific Blue™ anti-human CD8, Biotin anti-human CD8, APC anti-human CD8, Alexa Fluor® 700 anti-human CD8, Purified anti-human CD8 (Maxpar® Ready), Brilliant Violet 510™ anti-human CD8, Brilliant Violet 711™ anti-human CD8, Brilliant Violet 785™ anti-human CD8, Brilliant Violet 605™ anti-human CD8, PE/Dazzle™ 594 anti-human CD8, APC/Fire™ 750 anti-human CD8, Brilliant Violet 421™ anti-human CD8, TotalSeq™-A0046 anti-human CD8, TotalSeq™-C0046 anti-human CD8, Brilliant Violet 750™ anti-human CD8, TotalSeq™-B0046 anti-human CD8, Spark Blue™ 550 anti-human CD8, APC/Fire™ 810 anti-human CD8, PE/Fire™ 640 anti-human CD8, PE/Fire™ 700 anti-human CD8, TotalSeq™-D0046 anti-human CD8, GMP APC anti-human CD8, PE/Cyanine5 anti-human CD8 Antibody, Spark NIR™ 685 anti-human CD8

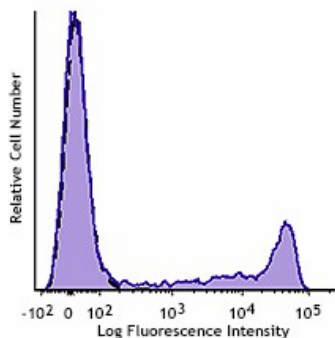

Typical results from human peripheral blood lymphocytes stained either with SK1 PE/Cyanine7 used at 5 µL/test (filled histogram) or with an isotype control (open histogram).

Symbols Glossary*

Symbol	Meaning	Symbol Title	Symbol No.	Symbol	Meaning	Symbol Title	Symbol No.
	Catalog number	Catalogue number	5.1.6		Indicates the need for the user to consult the instructions for use.	Consult instructions for use	5.4.3
	Indicates the temperature limits to which the medical device can be safely exposed.	Temperature limit	5.3.7		Indicates a medical device that needs protection from light sources.	Keep away from sunlight	5.3.2
	Indicates the upper limit of temperature to which the medical device can be safely exposed.	Upper limit of temperature	5.3.6		Indicates the date after which the medical device is not to be used.	Use-by date	5.1.4
	Indicates the medical device manufacturer.	Manufacturer	5.1.1		Indicates the authorized representative in the European Community.	Authorized representative in the European Community	5.1.2
	Indicates the manufacturer's batch code so that the batch or lot can be identified.	Batch code	5.1.5		Indicates a medical device that is intended to be used as an in vitro diagnostic medical device.	<i>In vitro</i> diagnostic medical device	5.5.1

* Symbol information is from EN ISO 15223-1:2016 Medical devices – Symbols to be used with medical device labels, labelling and information to be supplied – Part 1: General requirements

For research use only. Not for resale. BioLegend will not be held responsible for patent infringement or other violations that may occur with the use of our products.

*These products may be covered by one or more Limited Use Label Licenses (see the BioLegend Catalog or our website, www.biolegend.com/ordering#license). BioLegend products may not be transferred to third parties, resold, modified for resale, or used to manufacture commercial products, reverse engineer functionally similar materials, or to provide a service to third parties without written approval of BioLegend. By use of these products you accept the terms and conditions of all applicable Limited Use Label Licenses.

BioLegend Inc., 8999 BioLegend Way, San Diego, CA 92121 www.biolegend.com
 Toll-Free Phone: 1-877-Bio-Legend (246-5343) Phone: (858) 768-5800 Fax: (877) 455-9587