

## GMP APC/Fire™ 750 anti-human CD64 Antibody

<b>Catalog# / Size</b>	260368 / 100 tests
<b>Clone</b>	10.1
<b>Workshop</b>	VI MA36
<b>Other Names</b>	FcγRI, FcR I
<b>Isotype</b>	Mouse IgG1, κ
<b>Description</b>	CD64 is a 72 kD single chain type I glycoprotein also known as FcγRI and FcR I. CD64 is a member of the immunoglobulin superfamily and is expressed on monocytes/macrophages, dendritic cells, and activated granulocytes. The expression can be upregulated by IFN-γ stimulation. CD64 binds IgG immune complex. It plays a role in antigen capture, phagocytosis of IgG/antigen complexes, and antibody-dependent cellular cytotoxicity (ADCC).

### Product Details

---

<b>Reactivity</b>	Human
<b>Reported Reactivity</b>	Baboon, Capuchin Monkey, Chimpanzee, Squirrel Monkey
<b>Antibody Type</b>	Monoclonal
<b>Host Species</b>	Mouse
<b>Immunogen</b>	Human rheumatoid synovial fluid cells and fibronectin-purified monocytes.
<b>Formulation</b>	Phosphate-buffered solution, pH 7.2, containing True-Stain Monocyte Blocker™, 0.09% sodium azide and 0.2% (w/v) BSA (origin USA) and a stabilizer.
<b>Preparation</b>	The antibody was purified by affinity chromatography and conjugated with APC/Fire™ 750 under optimal conditions.
<b>Concentration</b>	200 µg/mL
<b>Storage &amp; Handling</b>	The antibody solution should be stored undiluted between 2°C and 8°C, and protected from prolonged exposure to light. <b>Do not freeze.</b>
<b>Application</b>	<a href="#">FC - Quality tested</a>
<b>Recommended Usage</b>	Each lot of this antibody is quality control tested by <a href="#">immunofluorescent staining with flow cytometric analysis</a> . For flow cytometric staining, the suggested use of this reagent is 5 µL per million cells in 100 µL staining volume or 5 µL per 100 µL of whole blood. It is recommended that the reagent be titrated for optimal performance for each application.  * APC/Fire™ 750 has a maximum excitation of 650 nm and a maximum emission of 774 nm.
<b>Excitation Laser</b>	Red Laser (633 nm)
<b>Application Notes</b>	Clone 10.1 recognizes the EC3 epitope of CD64. While both contain the EC3 domain, in-house testing suggests that clone 10.1 preferentially binds to CD64A (FcγRIA), but not CD64B (FcγRIB). Additional reported applications (for the relevant formats) include: blocking of human IgG3 and murine IgG2a binding to FcγRI <sup>2,5,6,11</sup> and immunohistochemical staining of acetone-fixed frozen tissue sections <sup>12</sup> .
<b>Application References</b>	
<b>(PubMed link indicates BioLegend citation)</b>	<ol style="list-style-type: none"><li>1. McMichael A, <i>et al.</i> Eds. 1987. Leucocyte Typing III. Oxford University Press. New York.</li><li>2. Schlossman S, <i>et al.</i> Eds. 1995. Leucocyte Typing V. Oxford University Press. New York. p. 874.</li><li>3. Kishimoto T, <i>et al.</i> Eds. 1997. Leucocyte Typing VI. Garland Publishing Inc. London.</li><li>4. Holl V, <i>et al.</i> 2004. <i>J. Immunol.</i> 173:6274.</li><li>5. Hober D, <i>et al.</i> 2002. <i>J. Gen. Virol.</i> 83:2169.</li><li>6. Cho HJ, <i>et al.</i> 2007. <i>Physiol Genomics</i> 149:60.</li></ol>

7. van Tits L, *et al.* 2005. *Arterioscler Thromb Vasc Biol.* 25:717. [PubMed](#)
8. Bruhns P, *et al.* 2008. *Blood* 113:3716. [PubMed](#)
9. Yoshino N, *et al.* 2000. *Exp. Anim. (Tokyo)* 49:97. (FC)
10. Carter DL, *et al.* 1999. *Cytometry* 37:41. (FC)
11. Dougherty GJ, *et al.* 1987. *Eur. J. Immunol.* 17:1453.
12. Blom AB, *et al.* 2003. *Arthritis Rheum.* 48(4):1002-14. (IHC)

#### Disclaimer

**GMP RUO Flow Cytometry Antibodies.** BioLegend GMP RUO fluorophore conjugated antibodies are manufactured in a dedicated GMP facility and compliant with ISO 13485:2016. For research use only. Not for use in diagnostic or therapeutic procedures. Our processes include:

- Batch-to-batch consistency
- Material traceability
- Documented procedures
- Documented employee training
- Equipment maintenance and monitoring records
- Lot-specific certificates of analysis
- Quality audits per ISO 13485:2016
- QA review of released products

## Antigen Details

---

<b>Structure</b>	Ig superfamily, type I glycoprotein, 72 kD
<b>Distribution</b>	Monocytes, macrophages, dendritic cells, activated granulocytes
<b>Function</b>	Phagocytosis, ADCC
<b>Ligand/Receptor</b>	IgG receptor
<b>Cell Type</b>	Dendritic cells, Granulocytes, Macrophages, Monocytes
<b>Biology Area</b>	Immunology, Innate Immunity
<b>Molecular Family</b>	CD Molecules, Fc Receptors
<b>Antigen References</b>	1. Hulett M, <i>et al.</i> 1994. <i>Adv. Immunol.</i> 57:1. 2. van de Winkel J, <i>et al.</i> 1993. <i>Immunol. Today</i> 14:215.
<b>Gene ID</b>	<a href="#">2209</a>

## Related Protocols

---

- [Cell Surface Flow Cytometry Staining Protocol](#)


## Other Formats

---

Biotin anti-human CD64, FITC anti-human CD64, PE anti-human CD64, Purified anti-human CD64, Alexa Fluor® 488 anti-human CD64, Alexa Fluor® 647 anti-human CD64, APC anti-human CD64, Pacific Blue™ anti-human CD64, Brilliant Violet 421™ anti-human CD64, PE/Cyanine7 anti-human CD64, PerCP/Cyanine5.5 anti-human CD64, APC/Cyanine7 anti-human CD64, Brilliant Violet 510™ anti-human CD64, Purified anti-human CD64 (Maxpar® Ready), PE/Dazzle™ 594 anti-human CD64, Brilliant Violet 605™ anti-human CD64, APC/Fire™ 750 anti-human CD64, TotalSeq™-A0162 anti-human CD64, Brilliant Violet 711™ anti-human CD64, Alexa Fluor® 700 anti-human CD64, Brilliant Violet 785™ anti-human CD64, TotalSeq™-C0162 anti-human CD64, Ultra-LEAF™ Purified anti-human CD64, TotalSeq™-B0162 anti-human CD64, TotalSeq™-D0162 anti-human CD64, GMP PE anti-human CD64, GMP FITC anti-human CD64, Brilliant Violet 650™ anti-human CD64, GMP PE/Dazzle™ 594 anti-human CD64

## Product Data

---

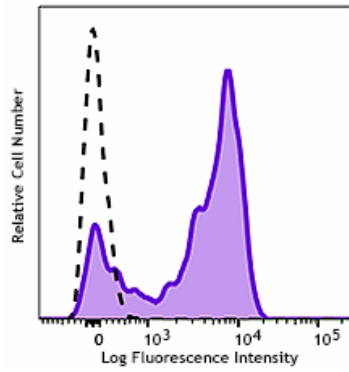

Typical results from human peripheral blood monocytes stained either with 10.1 APC/Fire™ 750 used at 5  $\mu$ L/test (filled histogram) or with an isotype control (open histogram).

For Research Use Only. Suitable for *ex vivo* cell processing. Not for injection or diagnostic or therapeutic use.

This product is supplied subject to the terms and conditions, including the limited license, located at [www.biolegend.com/terms](http://www.biolegend.com/terms) ("Terms") and may be used only as provided in the Terms. Without limiting the foregoing, BioLegend products may not be used for any Commercial Purpose as defined in the Terms, resold in any form, used in manufacturing, or reverse engineered, sequenced, or otherwise studied or used to learn its design or composition without express written approval of BioLegend. Regardless of the information given in this document, user is solely responsible for determining any license requirements necessary for user's intended use and assumes all risk and liability arising from use of the product. BioLegend is not responsible for patent infringement or any other risks or liabilities whatsoever resulting from the use of its products.

BioLegend, the BioLegend logo, and all other trademarks are property of BioLegend, Inc. or their respective owners, and all rights are reserved.

8999 BioLegend Way, San Diego, CA 92121 [www.biolegend.com](http://www.biolegend.com)  
Toll-Free Phone: 1-877-Bio-Legend (246-5343) Phone: (858) 768-5800 Fax: (877) 455-9587