

# BioLegend®

## LEGENDplex™ Multi-Analyte Flow Assay Kit

**Cat. No. 740891, Human Thrombosis Panel (10-plex) with Filter Plate**

Cat. No. 740892, Human Thrombosis Panel (10-plex) with V-Bottom Plate

**Cat. No. 740905, Human Thrombosis Panel (7-plex) with Filter Plate**

Cat. No. 740906, Human Thrombosis Panel (7-plex) with V-Bottom Plate

**Cat. No. 740907, Human Thrombosis Panel (3-plex) with Filter Plate**

Cat. No. 740908, Human Thrombosis Panel (3-plex) with V-Bottom Plate

**Please read the entire manual before running the assay.**

BioLegend.com

***It is highly recommended that this manual be read in its entirety before using this product. Do not use this kit beyond the expiration date.***

*For Research Purposes Only. Not for use in diagnostic or therapeutic procedures. Purchase does not include or carry the right to resell or transfer this product either as a stand-alone product or as a component of another product. Any use of this product other than the permitted use without the express written authorization of BioLegend is strictly prohibited.*

<b>Table of Contents</b>	<b>Page</b>
<b>Chapter 1: KIT DESCRIPTION</b> .....	3
Introduction.....	3
Principle of the Assay.....	4
Beads Usage.....	4
Storage Information.....	6
Materials Supplied.....	6
Materials to be Provided by the End-User.....	8
Precautions.....	9
<b>Chapter 2: ASSAY PREPARATION</b> .....	10
Sample Collection and Handling.....	10
Reagent Preparation.....	10
Standard Preparation.....	11
Sample Dilution.....	12
<b>Chapter 3: ASSAY PROCEDURE</b> .....	13
Performing the Assay Using a Filter Plate.....	13
Performing the Assay Using a V- bottom Plate.....	17
<b>Chapter 4: FLOW CYTOMETER SETUP</b> .....	20
<b>Chapter 5: DATA ACQUISITION AND ANALYSIS</b> .....	20
Data Acquisition.....	20
Data Analysis.....	21
<b>Chapter 6: ASSAY CHARACTERIZATION</b> .....	22
Representative Standard Curve.....	22
Assay Sensitivity.....	22
Cross-Reactivity.....	23
Accuracy.....	24
Linearity of Dilution.....	25

**LEGENDplex™ Human Thrombosis Panel**

Intra-Assay Precision..... 26  
Inter-Assay Precision..... 27  
Biological Samples..... 28

**TROUBLESHOOTING..... 31**

**PLATE MAP..... 37**

## Chapter 1: KIT DESCRIPTION

### Introduction

Thrombosis is the local coagulation of blood in response to vessel injury. Under normal conditions, platelets and fibrin assemble at the site of blood vessel injury and initiate clot formation to prevent further hemorrhage and re-establish hemostatic control. However, as the most common complication of the three major cardiovascular disorders, thrombosis also contributes to the pathologies of ischemic heart disease (acute coronary syndrome), stroke, and venous thromboembolism. Anti-coagulant therapy is often administered to patients considered at risk for thrombotic disorders in order to combat the progression of these diseases.

Amongst the several proteins involved in thrombosis, P-selectin, D-Dimer, tPA, PAI-1, Tissue Factor, and Factor IX act on fibrin at the site of injury, while IL-6, P-selectin, PSGL-1, sCD40L, and CXCL8 (IL-8) aid in the activation and assembly of platelets near the clot.

The Human Thrombosis Panel is a multiplex bead-based assay panel, using fluorescence–encoded beads suitable for use on various flow cytometers. This panel allows simultaneous quantification of 10 human proteins involved in thrombosis, including IL-6, P-selectin, D-Dimer, PSGL-1, tPA, sCD40L, PAI-1, Tissue Factor, Factor IX, and CXCL8 (IL-8). This assay panel provides broader dynamic ranges than traditional ELISA methods. The panel has been validated for use with tissue culture supernatant, serum, and plasma samples.

The LEGENDplex™ Human Thrombosis Panel is configured as shown below depending on sample types and required dilutions:

Catalog No.	Plex Size	Targets	Recommended Sample Type	Recommended Dilution Factor
740891 740892	10-plex	IL-6, P-selectin, D-Dimer, PSGL-1, tPA, sCD40L, PAI-1, Tissue Factor, Factor IX, CXCL8 (IL-8)	Tissue Culture	Varies
740905 740906	7-plex	P-selectin, D-Dimer, PSGL-1, tPA, sCD40L, PAI-1, Factor IX	Serum, Plasma	50
740907 740908	3-plex	IL-6, Tissue Factor, CXCL8 (IL-8)	Serum, Plasma	2

The LEGENDplex™ Human Thrombosis Panel is designed to allow flexible customization within the panel. Please visit [www.biolegend.com/legendplex](http://www.biolegend.com/legendplex) for more information on how to mix and match within the panel.

This assay is for research use only.

# LEGENDplex™ Human Thrombosis Panel

## Principle of the Assay

BioLegend's LEGENDplex™ assays are bead-based immunoassays using the same basic principle as sandwich immunoassays.

Beads are differentiated by size and internal fluorescence intensities. Each bead set is conjugated with a specific antibody on its surface and serves as the capture beads for that particular analyte. When a selected panel of capture beads is mixed and incubated with a sample containing target analytes, each analyte will bind to its specific capture beads. After washing, a biotinylated detection antibody cocktail is added, and each detection antibody in the cocktail will bind to its specific analyte bound on the capture beads, thus forming capture bead-analyte-detection antibody sandwiches. Streptavidin-phycoerythrin (SA-PE) is subsequently added, which will bind to the biotinylated detection antibodies, providing fluorescent signal intensities in proportion to the amount of bound analytes.


Since the beads are differentiated by size and internal fluorescence intensity on a flow cytometer, analyte-specific populations can be segregated and PE fluorescent signal quantified. The concentration of a particular analyte is determined using a standard curve generated in the same assay.

## Beads Usage

The Human Thrombosis Panel uses two sets of beads. Each set has a unique size that can be identified based on their forward scatter (FSC) and side scatter (SSC) profiles (Beads A and Beads B, Figure 1). Each bead set can be further resolved based on their internal fluorescence intensities. The internal dye can be detected using FL3, FL4, or APC channel, depending on the type of flow cytometer used. The smaller Beads A consists of 5 bead populations (A4, A5, A6, A7, A10) and the larger Beads B consists of 5 bead populations (B2, B4, B5, B6, B7) (Figure 2-3).

Using a total of 10 bead populations distinguished by size and internal fluorescent dye, the Human Thrombosis Panel allows simultaneous detection of 10 proteins in a single sample. Each analyte is associated with a particular bead set as indicated (Figures 2-3 and Table 1).

**Figure 1. Beads Differentiated by Size**

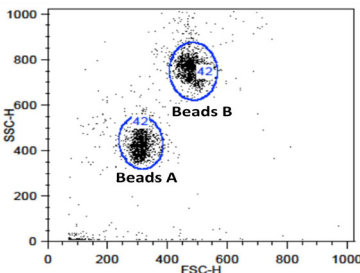

Beads A = smaller beads

Beads B = larger beads

Figure 2. Beads A Classification by FL4


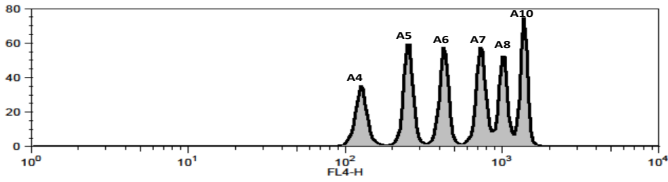


Figure 3. Beads B Classification by FL4

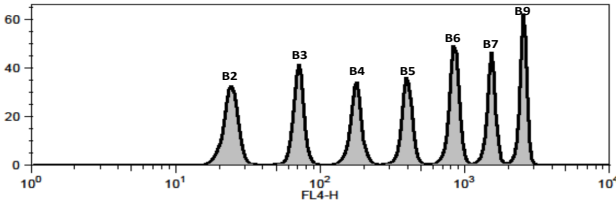

For Beads usage in the full panel as well as sample specific subpanels, please refer to Table 1 below.

Table 1. Panel Targets and Bead ID\*

Target	Bead ID	Thrombosis Panel (10-plex)	Thrombosis-Panel (7-plex)	Thrombosis Panel (3-plex)	Top Standard Concentrations
		Cat. # 740891 or 740892	Cat. # 740905 or 740906	Cat. # 740907 or 740908	
IL-6	A4	√		√	Note: The top standard concentrations of analytes in this panel were set at various concentrations, but may be subject to change from lot to lot (please visit <a href="http://biolegend.com/en-us/legendplex">biolegend.com/en-us/legendplex</a> to download a lot-specific certificate of analysis).
P-selectin	A5	√	√		
D-Dimer	A6	√	√		
PSGL-1	A7	√	√		
tPA	A10	√	√		
sCD40L	B2	√	√		
PAI-1	B4	√	√		
Tissue Factor	B5	√		√	
Factor IX	B6	√	√		
CXCL8 (IL-8)	B7	√		√	

Bead ID is used to associate a bead population to a particular analyte when using the LEGENDplex™ data analysis software program. For further information regarding the use of the program please visit [biolegend.com/en-us/legendplex](http://biolegend.com/en-us/legendplex).

# LEGENDplex™ Human Thrombosis Panel

## Storage Information

Recommended storage for all original kit components is between 2°C and 8°C. DO NOT FREEZE Pre-mixed Beads, Detection Antibodies or SA-PE.

- Once the standards have been sufficiently reconstituted, immediately transfer contents into polypropylene vials. DO NOT STORE RECONSTITUTED STANDARDS IN GLASS VIALS.
- Upon reconstitution, leftover top standard should be stored at  $\leq -70^{\circ}\text{C}$  for use within one month. Avoid multiple (>2) freeze-thaw cycles. Discard any leftover diluted standards.

## Materials Supplied

The LEGENDplex™ kit contains reagents for 100 tests, listed in the table below. When assayed in duplicate, this is enough for an 8-point standard curve and 40 samples.

Kit Components	Quantity	Volume	Part #
Setup Beads: PE Beads	1 vial	1 mL	77842
Setup Beads: Raw Beads	1 vial	1.8 mL	77844
<b>Capture Beads*</b> (see tables below for more information)	<b>Varies</b>	<b>Varies</b>	<b>Varies*</b>
<b>Human Thrombosis Panel Detection Antibodies</b>	<b>1 bottle</b>	<b>3.3 mL</b>	<b>Varies*</b>
<b>Human Thrombosis Panel Standard</b>	<b>1 vial</b>	<b>Lyophi- lized</b>	<b>Varies*</b>
LEGENDplex™ SA-PE	1 bottle	3.3 mL	77743
LEGENDplex™ Assay Buffer	1 bottle	25 mL	77562
Lyophilized Standard Reconstitution Buffer	1 vial	1 mL	75241
LEGENDplex™ Matrix B	1 vial	Lyophi- lized	77549
LEGENDplex™ Wash Buffer, 20X	1 bottle	25 mL	77564
Filter Plate** or V-bottom Plate***	1 plate		76187** or 76883***
Plate Sealers	4 sheets		78101

\* For full panel, premixed beads are provided ready-to-use. For 7-plex and 3-plex, individual beads are provided at 13X concentration. For Standard and Detection Antibodies, 10-plex uses part numbers, and 7-plex and 3-plex use catalog numbers (See tables below for details).

\*\* For kit with filter plate. \*\*\* For kit with V-bottom plate. Only one plate is provided for each kit.



## LEGENDplex™ Human Thrombosis Panel

### For Human Thrombosis Panel (10-plex):

Kit Components	Quantity	Volume	Part #
Human Thrombosis Panel Premixed Beads	1 bottle	3.3 mL	750000375
Human Thrombosis Panel Detection Antibodies	1 bottle	3.3 mL	750000373
Human Thrombosis Panel Standard	1 vial	Lyophi- lized	750000377

### For Human Thrombosis Panel (7-plex):

Kit Components	Quantity	Volume	Cat.#
LEGENDplex™ Human P-selectin Capture Bead A5, 13X	1 vial	270 µL	740896
LEGENDplex™ Human D-Dimer Capture Bead A6, 13X	1 vial	270 µL	740897
LEGENDplex™ Human PSGL-1 Capture Bead A7, 13X	1 vial	270 µL	741239
LEGENDplex™ Human tPA Capture Bead A10, 13X	1 vial	270 µL	740899
LEGENDplex™ Human sCD40L Capture Bead B2, 13X	1 vial	270 µL	740900
LEGENDplex™ Human PAI-1 Capture Bead B4, 13X	1 vial	270 µL	740901
LEGENDplex™ Human Factor IX Capture Bead B6, 13X	1 vial	270 µL	740903
LEGENDplex™ Human Thrombosis Panel Detection Antibodies	1 bottle	3.3 mL	740895
LEGENDplex™ Human Thrombosis Panel Standard	1 vial	Lyophi- lized	740894

### For Human Thrombosis Panel (3-plex):

Kit Components	Quantity	Volume	Cat.#
LEGENDplex™ Human IL-6 Capture Bead A4, 13X	1 vial	270 µL	740597
LEGENDplex™ Human Tissue Factor Capture Bead B5, 13X	1 vial	270 µL	740902
LEGENDplex™ Human CXCL8 (IL-8) Capture Bead B7, 13X	1 vial	270 µL	741240
LEGENDplex™ Human Thrombosis Panel Detection Antibodies	1 bottle	3.3 mL	740895
LEGENDplex™ Human Thrombosis Panel Standard	1 vial	Lyophi- lized	740894

Please refer to **Panel Targets and Bead ID (Table 1, page 5)**, to see which capture beads are included in each panel.

# LEGENDplex™ Human Thrombosis Panel

## Materials to be Provided by the End-User

- A flow cytometer equipped with two lasers (e.g., a 488 nm blue laser or 532 nm green laser and a 633-635 nm red laser) capable of distinguishing 575 nm and 660 nm or a flow cytometer equipped with one laser (e.g., 488 nm blue laser) capable of distinguishing 575 nm and 670 nm.

### Partial list of compatible flow cytometers:

Flow Cytometer	Reporter Channel	Reporter Emission	Classification Channel	Channel Emission	Compensation needed?
BD Accuri™ C6™	FL2	585 nm	FL4	675 nm	No*
BD FACSCanto, BD FACSCanto™II	PE	575 nm	APC	660 nm	No*
BD™ LSR, LSR II BD LSRFortessa™	PE	575 nm	APC	660 nm	No*
Gallios™	PE	575 nm	APC	660 nm	No*
CytoFLEX	PE	585 nm	APC	660 nm	No*
NovoCyte	PE	572 nm	APC	660 nm	No*
Attune™ NxT	PE	574 nm	APC	670 nm	No*

**\*Compensation is not required for the specified flow cytometers when set up properly.**

For setting up various flow cytometers, please visit: [www.biolegend.com/instrumentsetup](http://www.biolegend.com/instrumentsetup) and click on the **Instrument Setup** tab.

- Multichannel pipettes capable of dispensing 5 µL to 200 µL
- Reagent reservoirs for multichannel pipette
- Polypropylene microfuge tubes (1.5 mL)
- Laboratory vortex mixer
- Sonicator bath (e.g., Branson Ultrasonic Cleaner model #B200, or equivalent)
- Aluminum foil
- Absorbent pads or paper towels
- Plate shaker (e.g., Lab-Line Instruments model #4625, or equivalent)
- Tabletop centrifuges (e.g., Eppendorf centrifuge 5415 C, or equivalent)

## **If the assay is performed in a filter plate:**

- A vacuum filtration unit (Millipore MultiScreen® HTS Vacuum Manifold, cat # MSVMHTS00 or equivalent). Instructions on how to use the vacuum manifold can be found at the supplier's website.
- A vacuum source (mini vacuum pump or line vacuum, e.g., Millipore Vacuum Pump, catalog # WP6111560, or equivalent)
- If needed, additional Filter plates can be ordered from BioLegend (Cat# 740377 or 740378).

## **If the assay is performed in a V-bottom plate:**

- Centrifuge with a swinging bucket adaptor for microtiter plates (e.g., Beckman Coulter Allegra™ 6R Centrifuge with MICROPLUS CARRIER adaptor for GH3.8 and JS4.3 Rotors) .
- If needed, additional V-bottom plates can be ordered from BioLegend (Cat# 740379).

## **Precautions**

- All blood components and biological materials should be handled as potentially hazardous. Follow universal precautions as established by the Center for Disease Control and Prevention and by the Occupational Safety and Health Administration when handling and disposing of infectious agents.
- Sodium azide has been added to some reagents as a preservative. Although the concentrations are low, sodium azide may react with lead and copper plumbing to form highly explosive metal azides. On disposal, flush with a large volume of water to prevent azide build-up.
- Do not mix or substitute reagents from different kits or lots. Reagents from different manufacturers should not be used with this kit.
- Do not use this kit beyond its expiration date.
- SA-PE and beads are light-sensitive. Minimize light exposure.

## Chapter 2: ASSAY PREPARATION

### Sample Collection and Handling

#### **Preparation of Serum Samples:**

- Allow the blood to clot for at least 30 minutes and centrifuge for 20 minutes at 1,000 x *g*.
- Remove serum and assay immediately or aliquot and store samples at  $\leq -20^{\circ}\text{C}$ . Avoid multiple (>2) freeze/thaw cycles.
- When using frozen samples, it is recommended that samples be thawed completely, mixed and centrifuged to remove particulates prior to use.

#### **Preparation of Plasma Samples:**

- Plasma collection should be collected using an anti-coagulant (e.g., EDTA, Heparin, Citrate). Centrifuge for 20 minutes at 1,000 x *g* within 30 minutes of blood collection.
- Remove plasma and assay immediately, or aliquot and store samples at  $\leq -20^{\circ}\text{C}$ . Avoid multiple (>2) freeze/thaw cycles.
- When using frozen samples, it is recommended that samples be thawed completely, mixed well and centrifuged to remove particulates.

#### **Preparation of Tissue Culture Supernatant:**

- Centrifuge the sample to remove debris and assay immediately. If not possible, aliquot and store samples at  $\leq -20^{\circ}\text{C}$ . Avoid multiple (>2) freeze/thaw cycles.

### Reagent Preparation

#### **Preparation of Antibody-Immobilized Beads**

- **If pre-mixed beads are provided in the kit:**

Sonicate pre-mixed Beads bottle for 1 minute in a sonicator bath and then vortex for 30 seconds prior to use. If no sonicator bath is available, increase the vortexing time to 1 minute to completely resuspend the beads.

- **If individual beads (13X) are provided in the kit:**

The individual beads (13X) need to be combined with one another and diluted with Assay Buffer to create a 1X working solution of beads prior to use.

## LEGENDplex™ Human Thrombosis Panel

1. Sonicate each bead vial for 1 minute in a sonicator bath and then vortex for 30 seconds to completely resuspend the beads.
2. Calculate and prepare a 1X beads working solution based on the desired number of reactions and plex-size of your assay (i.e. the number of individual bead vials) following the steps described below.
  - A. Total volume ( $\mu\text{L}$ ) =  $30 \times$  (number of reactions)
  - B. Volume needed from each 13X beads vial ( $\mu\text{L}$ ) =  $2.3 \times$  (number of reactions)
  - C. Assay Buffer needed ( $\mu\text{L}$ ) =  $A - B \times$  (number of individual beads vials to be mixed)

Note: calculations for total volume include a 20% excess to account for any loss during pipetting.

### Example: to prepare 50 reactions for a 5-plex assay

- A. Total volume ( $\mu\text{L}$ ) =  $30 \times 50 = 1500 \mu\text{L}$
- B. Volume per beads vial needed ( $\mu\text{L}$ ) =  $2.3 \times 50 = 115 \mu\text{L}$
- C. Assay Buffer needed ( $\mu\text{L}$ ) =  $A - B \times$  (number of individual beads vials) =  $1500 - (115 \times 5) = 925 \mu\text{L}$

Combine 115  $\mu\text{L}$  of each beads vial (5 vials) with 925  $\mu\text{L}$  of assay buffer to get the desired final volume of 1500  $\mu\text{L}$  of 1X working solution of beads.

### Preparation of Wash Buffer

- Bring the 20X Wash Buffer to room temperature and mix to bring all salts into solution.
- Dilute 25 mL of 20X Wash Buffer with 475 mL deionized water. Store unused portions between 2°C and 8°C for up to one month.

### Preparation of Matrix B (for Serum or Plasma Samples Only)

- Add 5.0 mL LEGENDplex™ Assay Buffer to the bottle containing lyophilized Matrix B. Allow at least 15 minutes for complete reconstitution. Vortex to mix well. Leftover reconstituted Matrix B should be stored at  $\leq -70^\circ\text{C}$  for up to one month.

### Standard Preparation

1. Prior to use, reconstitute the lyophilized Human Thrombosis Standard with 250  $\mu\text{L}$  Lyophilized Standard Reconstitution Buffer.
2. Mix and allow the vial to sit at room temperature for 10 minutes, and then transfer the standard to an appropriately labeled polypropylene microcentrifuge tube. This will be used as the top standard C7.

## LEGENDplex™ Human Thrombosis Panel

**Note: The top standard concentrations of analytes in this panel were set at various concentrations, but may be subject to change from lot to lot (please visit [biolegend.com/en-us/legendplex](http://biolegend.com/en-us/legendplex) to download a lot-specific certificate of analysis).**

3. Label 6 polypropylene microcentrifuge tubes as C6, C5, C4, C3, C2 and C1, respectively.
4. Add 75  $\mu\text{L}$  of Assay Buffer to each of the six tubes. Prepare 1:4 dilution of the top standard by transferring 25  $\mu\text{L}$  of the top standard C7 to the C6 tube and mix well. This will be the C6 standard.
5. In the same manner, perform serial 1:4 dilutions to obtain C5, C4, C3, C2 and C1 standards (**see the table below using the top standard at 10,000  $\text{pg/mL}$  as an example**). Assay Buffer will be used as the 0  $\text{pg/mL}$  standard (C0).

Tube/Standard ID	Serial Dilution	Assay Buffer to add ( $\mu\text{L}$ )	Standard to add	Final Conc. ( $\text{pg/mL}$ )
C7	--	--	--	10,000
C6	1:4	75	25 $\mu\text{L}$ of C7	2,500
C5	1:16	75	25 $\mu\text{L}$ of C6	625
C4	1:64	75	25 $\mu\text{L}$ of C5	156.25
C3	1:256	75	25 $\mu\text{L}$ of C4	39.01
C2	1:1024	75	25 $\mu\text{L}$ of C3	9.77
C1	1:4096	75	25 $\mu\text{L}$ of C2	2.44
C0	--	75	--	0

### Sample Dilution

- For tissue culture supernatant samples, the levels of analyte can vary greatly from sample to sample. To test tissue culture supernatant samples, a preliminary experiment may be required to determine the appropriate dilution factor. If further dilution is desired, dilution should be done with corresponding fresh cell culture medium or Assay Buffer as a diluent to ensure accurate measurement.

For serum and plasma samples, follow panel specific dilution recommendations below:

- **For Human Thrombosis Panel (10-plex):**

Tissue culture samples are the recommended sample type. Follow the tissue culture supernatant dilution guidelines mentioned above.

- **For Human Thrombosis Panel (7-plex):**

Serum or plasma samples must be diluted 50-fold with Assay Buffer as described in the table below.

Sample	Dilution (1:50)	Final Dilution Fold
Serum, Plasma	2 $\mu$ L + 98 $\mu$ L (Assay Buffer)	50

If further dilution is desired, dilution should be done with Assay Buffer to ensure accurate measurement.

- **For Human Thrombosis Panel (3-plex):**

Serum or plasma samples must be diluted 2-fold with Assay Buffer as described in the table below.

Sample	Dilution (1:2)	Final Dilution Fold
Serum, Plasma	50 $\mu$ L + 50 $\mu$ L (Assay Buffer)	2

If further dilution is desired, dilution should be done with Matrix B to ensure accurate measurement.

## Chapter 3: ASSAY PROCEDURE

The LEGENDplex™ assay can be performed in a filter plate, or in a V-bottom plate.

- The Filter plate assay procedure requires a vacuum filtration unit for washing (see **Materials to be Provided by the End-User, page 8**). If you have performed bead-based multiplex assays before, your lab may already have the vacuum filtration unit set up.
- If the Filter plate assay procedure is not possible or if you prefer, the assay can be performed in a V-bottom plate.

### Performing the Assay Using a Filter Plate

- Allow all reagents to warm to room temperature (20-25°C) before use.
- Set the filter plate on an inverted plate cover at all times during assay setup and incubation steps, so that the bottom of the plate does not touch any surface. Touching a surface may cause leakage.
- Keep the plate upright during the entire assay procedure, including the washing steps, to avoid losing beads.

## LEGENDplex™ Human Thrombosis Panel

- The plate should be placed in the dark or wrapped with aluminum foil for all incubation steps.
- **Standards and samples should be run in duplicate and arranged on the plate in a vertical configuration convenient for data acquisition and analysis (as shown in attached PLATE MAP, page 37). Be sure to load standards in the first two columns. If an automation device is used for reading, the orientation and reading sequence should be carefully planned.**

1. Pre-wet the plate by adding 100 µL of LEGENDplex™ 1X Wash Buffer to each well and let it sit for 1 minute at room temperature. To remove the excess volume, place the plate on the vacuum manifold and apply vacuum. Do not exceed 10" Hg of vacuum. Vacuum until wells are drained (5-10 seconds). Blot excess Wash Buffer from the bottom of the plate by pressing the plate on a stack of clean paper towels. Place the plate on top of the inverted plate cover.
2. **For measuring tissue culture supernatant samples in 10-plex, and serum or plasma samples in 7-plex,** load the plate as shown in the table below (in the order from left to right):

	Cell Culture Medium or Assay Buffer	Standard	Sample*
Standard Wells	25 µL	25 µL	---
Sample wells	25 µL	---	25 µL

**For measuring serum or plasma samples in 3-plex,** load the plate as shown in the table below (in the order from left to right):

	Matrix B	Assay Buffer	Standard	Sample*
Standard Wells	25 µL	---	25 µL	---
Sample wells	---	25 µL	---	25 µL

\*See **Sample Dilution on page 12**

3. Vortex mixed beads bottle for 30 seconds. Add 25 µL of mixed beads to each well. The volume should be 75 µL in each well after beads addition. (Note: During addition of the beads, shake mixed beads bottle intermittently to avoid bead settling).
4. Seal the plate with a plate sealer. **To avoid plate leaking, do not apply positive pressure to the sealer when sealing the plate.** Wrap the entire plate, including the inverted plate cover, with aluminum foil. Place the plate on a plate shaker, secure it with a rubber band and shake at approximate 500 rpm for 2 hours at room temperature.
5. **Do not invert the plate!** Place the plate on the vacuum manifold and apply vacuum as before in Step 1. Add 200 µL of 1X Wash Buffer to each well.



## LEGENDplex™ Human Thrombosis Panel

Remove Wash Buffer by vacuum filtration. Blot excess Wash Buffer from the bottom of the plate with an absorbent pad or paper towels. Repeat this washing step once more.

6. Add 25  $\mu$ L of Detection Antibodies to each well.
7. Seal the plate with a fresh plate sealer. Wrap the entire plate, including the inverted plate cover, with aluminum foil. Place the plate on a plate shaker and shake at approximately 500 rpm for 1 hour at room temperature.
8. **Do not vacuum!** Add 25  $\mu$ L of SA-PE to each well directly.
9. Seal the plate with a fresh plate sealer. Wrap the entire plate, including the inverted plate cover, with aluminum foil. Place the plate on a plate shaker and shake at approximate 500 rpm for 30 minutes at room temperature.
10. Repeat step 5 above.
11. Add 150  $\mu$ L of 1X Wash Buffer to each well. Resuspend the beads on a plate shaker for 1 minute.
12. Read samples on a flow cytometer, preferably within the same day of the assay (Note: Prolonged sample storage can lead to reduced signal).


If the flow cytometer is equipped with an autosampler, read the plate directly using the autosampler. **Please be sure to program the autosampler to resuspend beads in the well immediately before taking samples. The probe height may need to be adjusted when using an autosampler.**

If an autosampler is not available, the samples can be transferred from the filter plate to micro FACS (or FACS) tubes and read manually.

# LEGENDplex™ Human Thrombosis Panel

## Assay Procedure Summary for Filter Plate

Add 100  $\mu$ L 1X Wash Buffer to filter plate wells  
Vacuum to remove excess buffer

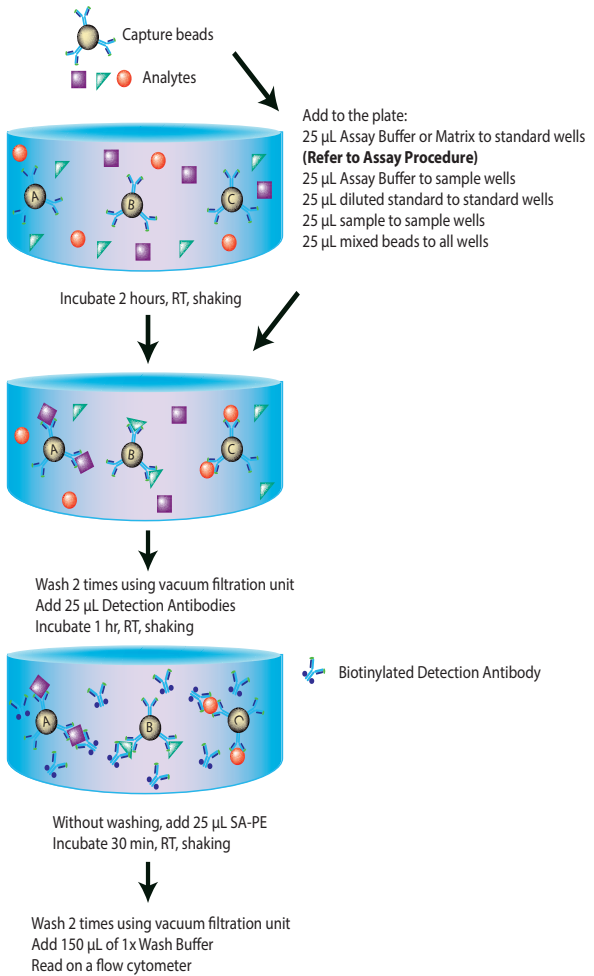

**Performing the Assay Using a V-bottom Plate**

- Allow all reagents to warm to room temperature (20-25°C) before use.
- Keep the plate upright during the entire assay procedure, including the washing steps, to avoid losing beads.
- The plate should be placed in the dark or wrapped with aluminum foil for all incubation steps.
- Standards and samples should be run in duplicate and arranged on the plate in a vertical configuration convenient for data acquisition and analysis (as shown in attached PLATE MAP, page 37). Be sure to load standards in the first two columns. If an automation device is used for reading, the orientation and reading sequence should be carefully planned.

1. **For measuring tissue culture supernatant samples in 10-plex, and serum or plasma samples in 7-plex**, load the plate as shown in the table below (in the order from left to right):

	Cell Culture Medium or Assay Buffer	Standard	Sample*
Standard Wells	25 µL	25 µL	---
Sample wells	25 µL	---	25 µL

**For measuring serum or plasma samples in 3-plex**, load the plate as shown in the table below (in the order from left to right):

	Matrix B	Assay Buffer	Standard	Sample*
Standard Wells	25 µL	---	25 µL	---
Sample wells	---	25 µL	---	25 µL

\*See **Sample Dilution on page 12**

2. Vortex mixed beads for 30 seconds. Add 25 µL of mixed beads to each well. The total volume should be 75 µL in each well after beads addition. (Note: During beads addition, shake mixed beads bottle intermittently to avoid bead settling).
3. Seal the plate with a plate sealer. Cover the entire plate with aluminum foil to protect the plate from light. Shake at 800 rpm on a plate shaker for 2 hours at room temperature (**Depending on the shaker, the speed may need to be adjusted. The optimal speed is one that is high enough to keep beads in suspension during incubation, but not too high that it may cause sample to spill from the wells**).
4. Centrifuge the plate at 1050 rpm (~250 g) for 5 minutes, using a swinging bucket rotor (G.H 3.8) with microplate adaptor (Please refer to **Materials to be Provided by the End-User, page 8**). Do not use excessive centrifugation

## LEGENDplex™ Human Thrombosis Panel

speed as it may make it harder to resuspend beads in later steps.

**Make sure the timer of the centrifuge works properly and standby to make sure the centrifuge reaches preset speed.**

5. Immediately after centrifugation, dump the supernatant into a biohazard waste container by quickly inverting and flicking the plate **in one continuous and forceful motion**. The beads pellet may or may not be visible after dumping the supernatant. Loss of beads should not be a concern as the beads will stay in the tip of the well nicely. Blot the plate on a stack of clean paper towel and drain the remaining liquid from the well as much as possible. Be careful not to disturb the bead pellet.

Alternatively, removal of the supernatant may be completed using a multichannel pipette set at 75  $\mu\text{L}$ . Try to remove as much liquid as possible without removing any beads. Be sure to change pipette tips between each row or column.


6. Wash the plate by dispensing 200  $\mu\text{L}$  of 1X Wash Buffer into each well and incubate for one minute. Repeat step 4 and 5 above. A second wash is optional, but may help reduce background.
7. Add 25  $\mu\text{L}$  of Detection Antibodies to each well.
8. Seal the plate with a new plate sealer. Cover the entire plate with aluminum foil to protect the plate from light. Shake at 800 rpm on a plate shaker for 1 hour at room temperature.
9. **Do not wash the plate!** Add 25  $\mu\text{L}$  of SA-PE to each well directly.
10. Seal the plate with a new plate sealer. Wrap the entire plate with aluminum foil and shake the plate on a plate shaker at approximate 800 rpm for 30 minutes at room temperature.
11. Repeat step 4 and 5.
12. (This washing step is optional but helps to reduce the background.) Wash the plate by dispensing 200  $\mu\text{L}$  of 1X Wash Buffer into each well and incubate for one minute. Repeat step 4 and 5 above.
13. Add 150  $\mu\text{L}$  of 1X Wash Buffer to each well. Resuspend the beads by pipetting.
14. Read samples on a flow cytometer, preferably within the same day of the assay (Note: Prolonged sample storage can lead to reduced signal).

If the flow cytometer is equipped with an autosampler, the samples can be read directly. **Please be sure to program the autosampler to resuspend beads in the well immediately before taking samples. The probe height may need to be adjusted when using an autosampler.**

If an autosampler is not available, the samples can be transferred from the plate to micro FACS (or FACS) tubes and read manually.

# LEGENDplex™ Human Thrombosis Panel

## Assay Procedure Summary for V-bottom Plate

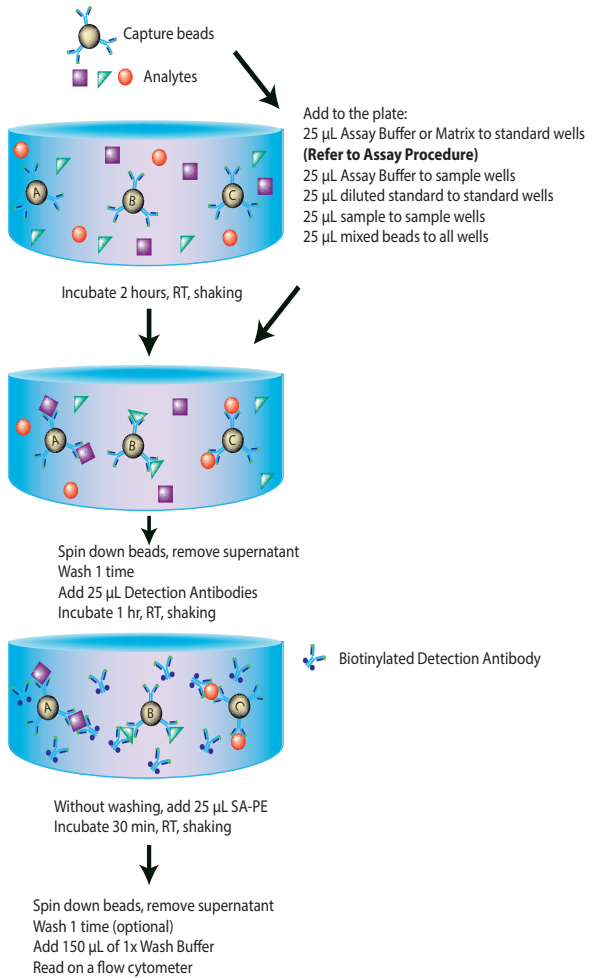

## Chapter 4: FLOW CYTOMETER SETUP

In order to generate reliable data, the flow cytometer must be set up properly before data acquisition.

The setup instructions have been removed from this manual and uploaded onto our website to save paper.

To access the setup instructions, please visit: [www.biolegend.com/legendplex](http://www.biolegend.com/legendplex) and click on the **Instrument Setup** tab.

## Chapter 5: DATA ACQUISITION AND ANALYSIS

### Data Acquisition

1. Before reading samples, make sure that the flow cytometer is set up properly.
2. Create a new template or open an existing template (for details on how to create a cytometer-specific template, please refer to the Flow Cytometer Setup Guide).
3. Vortex each sample for 5 seconds before analysis.
4. Set the flow rate to low. Set the number of beads to be acquired to about 300 per analyte (e.g., acquire 900 beads for a 3-plex assay or 3,000 beads for a 10-plex assay). Do not set to acquire total events as samples may contain large amounts of debris. Instead, create a large gate to include both Beads A and Beads B (gate A+B) and set to acquire the number of events in gate A + B. This will exclude majority of the debris.

Note: Do not acquire too few or too many beads. Too few beads acquired may result in high CVs and too many beads acquired may result in slow data analysis later.

5. Read samples.

When reading samples, set the flow cytometer to setup mode first and wait until bead population is stabilized before recording or switching to acquisition mode.

To simplify data analysis using the LEGENDplex™ Data Analysis Software, read samples in the same order as shown on the PLATE MAP attached at the end of the manual. For an in-plate assay, read column by column (A1, B1, C1...A2, B2, C2...).

## LEGENDplex™ Human Thrombosis Panel

When naming data files, try to use simple names with a consecutive numbering for easy data analysis (e.g. for standards, C0.001, C0.002, C1.003, C1.004, C2.005, C2.006, C3.007, C3.008, ... C7.015, C7.016; for samples, S1.017, S1.018, S2.019, S2.020, S3.021, S3.022...)

Store all FCS files in the same folder for each assay. If running multiple assays, create a separate folder for each assay.

6. Proceed to data analysis using LEGENDplex™ Data Analysis Software when data acquisition is completed.


### **Data Analysis**

- The assay FCS files should be analyzed using BioLegend's LEGENDplex™ data analysis software. The program is offered free of charge with the purchase of any LEGENDplex™ assay. For further information regarding access to, and use of the program please visit [biolegend.com/en-us/legendplex](http://biolegend.com/en-us/legendplex).

## Chapter 6: ASSAY CHARACTERIZATION

### Representative Standard Curve

This standard curve was generated using the LEGENDplex™ Human Thrombosis Panel for demonstration purposes only. A standard curve must be run with each assay.

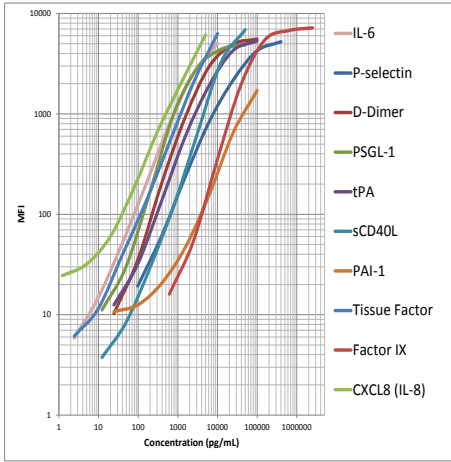

### Assay Sensitivity

The assay sensitivity is the theoretical limit of detection calculated using the LEGENDplex™ Data Analysis Software by applying a 5-parameter curve fitting algorithm. Assay Sensitivity presented here is  $\leq$  Mean LOD + 2x STDEV LOD.

Analyte	Sensitivity in Assay Buffer (pg/mL)	Sensitivity in Matrix (pg/mL)
IL-6	1.53	3.00
P-selectin	49.00	N/A
D-Dimer	12.00	N/A
PSGL-1	20.73	N/A
tPA	12.00	N/A
sCD40L	12.32	N/A
PAI-1	66.45	N/A
Tissue Factor	4.80	4.84



## LEGENDplex™ Human Thrombosis Panel

Factor IX	305.00	N/A
CXCL8 (IL-8)	6.06	6.00

### Cross-Reactivity

Target proteins were tested individually at the indicated concentrations below using the LEGENDplex™ Human Thrombosis Panel, with negligible cross-reactivity observed for non-intended targets.

Analyte	Conc. (ng/mL)	Analyte	Conc. (ng/mL)
IL-6	100	sCD40L	500
P-selectin	4,000	PAI-1	1,000
D-Dimer	1,000	Tissue Factor	100
PSGL-1	500	Factor IX	25,000
tPA	1,000	CXCL8 (IL-8)	50

The following recombinant proteins were tested individually at 50 ng/mL. No or negligible cross-reactivity was found.

Human				Mouse	
uPA	CXCL10	Myoglobin	ICAM-1	LIGHT	IL-6
Vitronectin	IFN-γ	MRP8/14	VCAM-1	TNF-α	L-selectin
Serpin E2	GM-CSF	NGAL	MMP-9	PIGF	PSGL-1
Serpin F2	IL-4	CRP	Cystatin C	IL-18	
L-selectin	Protein C	MMP-2	ST2	IL-10	
E-selectin	Protein S	OPN	RAGE	CCL2	
IL-12(p70)	Protein Z	MPO	TIE-2	PTX3	
IL-17A	Factor VII	SAA	TIE-1	VWF	
IL-1β	Factor X	IGFBP4	Flt-1	*Fibrinogen	

\*tested at 80 μg/mL with negligible cross-reactivity observed.

**Accuracy (Spike Recovery)**

For spike recovery in tissue culture supernatant, serum, and plasma, target proteins with known concentrations in the assay range were spiked into sample. The spiked samples were then assayed, and the measured concentrations were compared with the expected values.

Analyte	% of Spike Recovery		
	Tissue Culture Supernatant (N= 2)	Serum (N = 4)	Plasma (N = 4)
IL-6	86%	86%	74%
P-selectin	88%	99%	94%
D-Dimer	77%	90%	77%
PSGL-1	89%	80%	79%
tPA	86%	96%	84%
sCD40L	82%	104%	85%
PAI-1	81%	153%	107%
Tissue Factor	92%	66%	60%
Factor IX	89%	79%	68%
CXCL8 (IL-8)	90%	76%	69%

**Linearity of Dilution**

- Tissue culture samples were spiked with target proteins with known concentrations in the assay range, then serially diluted 1:2, 1:4, 1:8 with Assay Buffer and assayed.
- Serum and plasma samples were first diluted 50-fold with Assay Buffer, then serially diluted 1:2, 1:4, 1:8 with Assay Buffer and assayed with the LEGENDplex™ Human Thrombosis Panel (7-plex) kit.
- Serum and plasma samples were spiked\* with target proteins with known concentrations in the assay range, then serially diluted 1:2, 1:4, 1:8 with Matrix B and assayed with the LEGENDplex™ Human Thrombosis Panel (3-plex) kit.

The measured concentrations of serially diluted samples were then compared with the concentration of the lowest dilution based on serial dilution factor used.

Analyte	% Linearity		
	Tissue Culture Supernatant (N= 2)	Serum (N = 4)	Plasma (N = 4)
IL-6	99%	104%*	111%*
P-selectin	120%	85%	79%
D-Dimer	91%	128%	116%
PSGL-1	108%	117%	100%
tPA	112%	100%	92%
sCD40L	103%	104%	96%
PAI-1	100%	106%	141%
Tissue Factor	104%	118%*	118%*
Factor IX	103%	117%	105%
CXCL8 (IL-8)	94%	117%*	134%*

# LEGENDplex™ Human Thrombosis Panel

## Intra-Assay Precision

Two samples with different concentrations of each target protein were analyzed in one assay with 16 replicates per sample. The intra-assay precision is shown below.

Analyte	Sample	Mean (pg/mL)	STDEV	%CV
IL-6	Sample 1	39.19	4.79	12%
	Sample 2	159.35	15.17	10%
P-selectin	Sample 1	1465.24	148.32	10%
	Sample 2	5991.00	430.50	7%
D-Dimer	Sample 1	419.92	20.94	5%
	Sample 2	1461.80	57.22	4%
PSGL-1	Sample 1	176.57	11.48	7%
	Sample 2	747.12	34.80	5%
tPA	Sample 1	327.25	35.55	11%
	Sample 2	1282.59	96.94	8%
sCD40L	Sample 1	52.18	5.66	11%
	Sample 2	282.27	20.24	7%
PAI-1	Sample 1	236.32	33.72	14%
	Sample 2	1161.91	104.72	9%
Tissue Factor	Sample 1	40.85	2.77	7%
	Sample 2	165.48	10.91	7%
Factor IX	Sample 1	9824.51	600.62	6%
	Sample 2	35580.47	1719.79	5%
CXCL8 (IL-8)	Sample 1	22.63	3.88	17%
	Sample 2	88.49	11.11	13%

**Inter-Assay Precision**

Two samples with different concentrations of each target protein were analyzed in four independent assays with four replicates per sample. The inter-assay precision is shown below.

Analyte	Sample	Mean (pg/mL)	STDEV	%CV
IL-6	Sample 1	54.46	5.64	10%
	Sample 2	220.67	24.89	11%
P-selectin	Sample 1	1862.78	271.08	15%
	Sample 2	6970.12	743.18	11%
D-Dimer	Sample 1	488.73	43.63	9%
	Sample 2	1692.14	114.96	7%
PSGL-1	Sample 1	220.70	28.37	13%
	Sample 2	879.07	57.72	7%
tPA	Sample 1	414.82	61.06	15%
	Sample 2	1270.46	129.21	10%
sCD40L	Sample 1	85.19	14.09	17%
	Sample 2	214.84	26.78	12%
PAI-1	Sample 1	355.51	74.11	21%
	Sample 2	1380.65	145.08	11%
Tissue Factor	Sample 1	49.40	8.10	16%
	Sample 2	198.89	18.56	9%
Factor IX	Sample 1	12726.00	1793.84	14%
	Sample 2	45100.38	3134.20	7%
CXCL8 (IL-8)	Sample 1	28.38	6.25	22%
	Sample 2	120.82	18.57	15%

# LEGENDplex™ Human Thrombosis Panel

## Biological Samples

### Serum and Plasma (Samples are paired)

Normal human serum samples (N = 20) were tested for levels of the target proteins. The concentrations measured are shown below.

Analyte	Range (pg/mL)	% of Detectable	Median (pg/mL)
IL-6	ND-48.18	85%	7.70
P-selectin	4,637.00-17,391.50	100%	8506.25
D-Dimer	10,742.50-48,237.50	100%	27421.25
PSGL-1	5,594.00-10,109.50	100%	7474.00
tPA	1,914.00-8,699.00	100%	5243.75
sCD40L	1,346.00-19,975.50	100%	6701.25
PAI-1	102,847.00-335,181.00	100%	192170.75
Tissue Factor	ND-172.72	95%	3.49
Factor IX	689,121.00-2,336,085.00	100%	1263790.50
CXCL8 (IL-8)	ND-81.64	75%	4.32

ND = Non-detectable

Normal human plasma samples (n=20) were tested for levels of the target proteins. The concentrations measured are shown below.

Analyte	Range (pg/mL)	% of Detectable	Median (pg/mL)
IL-6	1.48-52.10	100%	13.52
P-selectin	1,930.00-216,632.50	100%	9921.25
D-Dimer	6,759.00-582,389.00	100%	26013.25
PSGL-1	3,574.00-12,050.00	100%	6509.25
tPA	1,588.00-8,392.00	100%	4647.50
sCD40L	2,338.00-76,602.50	100%	5459.25
PAI-1	86,070.50-659,364.00	100%	202673.75
Tissue Factor	5.94-200.90	100%	38.67
Factor IX	835,986.50-2,475,181.00	100%	1311631.50
CXCL8 (IL-8)	11.34-133.74	100%	43.23

**Deep Vein Thrombosis Patient Samples**

Deep Vein Thrombosis patient serum samples (N = 4) were purchased from a commercial source and tested for levels of the target proteins. The concentrations measured are shown below.

Analyte	Range (pg/mL)	% of Detectable	Median (pg/mL)
IL-6	1.98-19.46	100%	8.07
P-selectin	12,690.50-24,455.50	100%	16461.50
D-Dimer	57,275.00-2,073,943.00	100%	153606.00
PSGL-1	8,607.50-14,937.00	100%	11065.50
tPA	3,617.00-14,711.00	100%	10925.25
sCD40L	3,367.50-38,146.50	100%	7492.00
PAI-1	79,467.50-520,762.00	100%	261485.50
Tissue Factor	2.90-8.28	100%	4.92
Factor IX	1,032,462.00-2,546,994.50	100%	1932488.50
CXCL8 (IL-8)	ND-1.94	25%	ND

ND = Non-detectable

Deep Vein Thrombosis patient plasma samples (N = 3) were purchased from a commercial source and tested for levels of the target proteins. The concentrations measured are shown below.

Analyte	Range (pg/mL)	% of Detectable	Median (pg/mL)
IL-6	4.54-11.60	100%	11.28
P-selectin	23,035.00-119,257.00	100%	112617.00
D-Dimer	18,134.50-214,489.50	100%	77959.00
PSGL-1	6,474.00-10,096.00	100%	9487.50
tPA	6,918.50-14,584.00	100%	7137.50
sCD40L	6,801.00-31,647.50	100%	24404.50
PAI-1	270,732.50-502,254.00	100%	332720.50
Tissue Factor	5.34-18.66	100%	5.90
Factor IX	1,024,499.50-2,344,930.00	100%	1870466.00
CXCL8 (IL-8)	2.76-26.14	100%	8.96

# LEGENDplex™ Human Thrombosis Panel

## Cell Culture Supernatant

Human PBMCs ( $1 \times 10^6$  cells/mL) were cultured under various stimulation conditions (Poly I:C, 50  $\mu\text{g/mL}$ ; R848, 2  $\mu\text{g/mL}$ ), using unstimulated cells as a control. Supernatants were collected 9 hours after stimulation, and assayed with the LEGENDplex™ Human Thrombosis Panel (10-plex) kit. The results (in  $\text{pg/mL}$ ) are summarized below.

Analyte	PBMC Control	PBMC + Poly I:C	PBMC + R848
IL-6	376.22	4584.46	13,211.06
P-selectin	5507.72	5752.48	4145.33
D-Dimer	37,795.89	25,963.65	66,937.84
PSGL-1	243.06	70.95	290.60
tPA	ND	ND	43.41
sCD40L	352.24	473.12	276.51
PAI-1	21,448.42	23,282.55	28,182.21
Tissue Factor	5.03	ND	7.71
Factor IX	ND	ND	ND
CXCL8 (IL-8)	3690.90	512.95	5955.90

ND = Non-detectable

Human umbilical endothelial cells (HUVECs) were cultured under various stimulation conditions (LPS, 1  $\mu\text{g/mL}$ ; 50  $\text{ng/mL}$  IFN- $\gamma$  with 200  $\text{ng/mL}$  TNF- $\alpha$ ), using unstimulated cells as a control. Supernatants were collected 24 hours after stimulation, and assayed with the LEGENDplex™ Human Thrombosis Panel (10-plex) kit. The results (in  $\text{pg/mL}$ ) are summarized below.

Analyte	HUVEC Control	HUVEC + LPS	HUVEC + IFN- $\gamma$ + TNF- $\alpha$
IL-6	133.36	668.05	2117.28
P-selectin	ND	ND	ND
D-Dimer	ND	ND	ND
PSGL-1	56.84	112.28	126.33
tPA	37836.57	42661.37	38536.01
sCD40L	ND	6.72	11.71
PAI-1	>100,000	>100,000	>100,000
Tissue Factor	0.94	1.10	1.35
Factor IX	ND	ND	ND
CXCL8 (IL-8)	905.95	4166.79	6061.28

ND = Non-detectable



## TROUBLESHOOTING

Problem	Possible Cause	Solution
Bead population shifting upward or downward during acquisition	The strong PE signal from high concentration samples or standards may spill over to classification Channel (e.g., FL3/FL4/APC) and mess up the bead separation.	Optimize instrument settings using Kit Setup Beads, and make appropriate compensation between channels.
Filter plate will not vacuum or some wells clogged	Vacuum pressure is insufficient or vacuum manifold does not seal properly.	Increase vacuum pressure such that 0.2 mL buffer can be suctioned in 3-5 seconds. Clean the vacuum manifold and make sure no debris on the manifold. Press down the plate on the manifold to make a good seal.
	Samples have insoluble particles or sample is too viscous (e.g., serum and plasma samples)	<p>Centrifuge samples just prior to assay setup and use supernatant. If high lipid content is present, remove lipid layer after centrifugation. Sample may need dilution if too viscous.</p> <p>If some wells are still clogged during washing, try the following:</p> <ol style="list-style-type: none"> <li>1). Add buffer to all the wells, pipette up and down the clogged wells and vacuum again.</li> <li>2). Use a piece of clean wipe, wipe the under side of the clogged wells and vacuum again.</li> <li>3). Take a thin needle (e.g., insulin needle), while holding the plate upward, poke the little hole under each of the clogged wells and vacuum again. Do not poke too hard or too deep as it may damage the filter and cause leaking.</li> </ol>
	Filter plate was used without pre-wet.	Pre-wet plate with wash buffer before running the assay.

# LEGENDplex™ Human Thrombosis Panel

Insufficient bead count or slow reading	Beads inappropriately prepared	Sonicate bead vials and vortex just prior to addition. Agitate mixed beads intermittently in reservoir while pipetting this into the plate.
	Samples cause beads aggregation due to particulate matter or viscosity.	Centrifuge samples just prior to assay setup and use supernatant. If high lipid content is present, remove lipid layer after centrifugation. Sample may need dilution if too viscous.
	Beads were lost during washing for in-tube assay	Make sure beads are spun down by visually check the pellet (beads are in light blue or blue color). Be very careful when removing supernatant during washing.
	Probe might be partially clogged.	Sample probe may need to be cleaned, or if needed, probe should be removed and sonicated.
Plate leaked	Vacuum pressure set too high	Adjust vacuum pressure such that 0.2 mL buffer can be suctioned in 3-5 seconds. Do not exceed 10" Hg of vacuum.
	Plate set directly on table or absorbent towels during incubations or reagent additions	Set plate on plate holder or raised edge so bottom of filter is not touching any surface.
	Liquid present on the under side of the plate after vacuum	After washing, press down plate firmly on a stack of clean paper towels to dry the underside of the plate.
	Pipette touching and damaged plate filter during additions.	Pipette to the side of wells.
High Background	Background wells were contaminated	Avoid cross-well contamination by changing tips between pipetting when performing the assay using a multichannel pipette.
	Insufficient washes	The background may be due to non-specific binding of SA-PE. Increase number of washes.
Debris (FSC/SSC) during sample acquisition	Debris or platelet may exist in sample solution.	Centrifuge samples before analyzing samples. Remove platelet as much as possible.

## LEGENDplex™ Human Thrombosis Panel

Variation between duplicate samples	Beads aggregation	Sonicate and vortex the Beads prior to use.
	Multichannel pipette may not be calibrated or inconsistent pipetting	Calibrate Pipette. Ensure good pipetting practice. Prime pipette before use may help.
	Plate washing was not uniform	Make sure all reagents are vacuumed out completely in all wash steps.
	Samples may contain particulate matters.	Centrifuge samples just prior to assay setup and use supernatant. If high lipid content is present, remove lipid layer after centrifugation. Sample may need dilution if too viscous.
Low or poor standard curve signal	The standard was incorrectly reconstituted, stored or diluted	Follow the protocol to reconstitute, store and dilute standard. Double check your calculation.
	Wrong or short incubation time	Ensure the time of all incubations was appropriate.
Signals too high, standard curves saturated	PMT value for FL2/PE set too high	Make sure the PMT setting for the reporter channel is appropriate
	Plate incubation time was too long	Use shorter incubation time.
Sample readings are out of range	Samples contain no or below detectable levels of analyte	Make sure the experiment to generate the samples worked. Use proper positive controls.
	Samples concentrations higher than highest standard point.	Dilute samples and analyze again.
	Standard curve was saturated at higher end of curve.	Make sure the PMT setting for the reporter channel is appropriate. Use shorter incubation time if incubation time was too long
Missed beads populations during reading, or distribution is unequal	Sample may cause some beads to aggregate.	Centrifuge samples just prior to assay setup and use supernatant. If high lipid content is present, remove lipid layer after centrifugation. Sample may need dilution if too viscous.
	Beads populations are not mixed properly	Make sure all bead populations are mixed. and in similar numbers.







**PLATE MAP (for in-plate assay)**

	1	2	3	4	5	6	7	8	9	10	11	12
<b>A</b>	<b>C0</b>	<b>C4</b>	Sample1	Sample5	Sample 9	Sample 13	Sample 17	Sample 21	Sample 25	Sample 29	Sample 33	Sample 37
<b>B</b>	<b>C0</b>	<b>C4</b>	Sample1	Sample5	Sample 9	Sample 13	Sample 17	Sample 21	Sample 25	Sample 29	Sample 33	Sample 37
<b>C</b>	<b>C1</b>	<b>C5</b>	Sample2	Sample6	Sample 10	Sample 14	Sample 18	Sample 22	Sample 26	Sample 30	Sample 34	Sample 38
<b>D</b>	<b>C1</b>	<b>C5</b>	Sample2	Sample6	Sample 10	Sample 14	Sample 18	Sample 22	Sample 26	Sample 30	Sample 34	Sample 38
<b>E</b>	<b>C2</b>	<b>C6</b>	Sample3	Sample7	Sample 11	Sample 15	Sample 19	Sample 23	Sample 27	Sample 31	Sample 35	Sample 39
<b>F</b>	<b>C2</b>	<b>C6</b>	Sample3	Sample7	Sample 11	Sample 15	Sample 19	Sample 23	Sample 27	Sample 31	Sample 35	Sample 39
<b>G</b>	<b>C3</b>	<b>C7</b>	Sample4	Sample8	Sample 12	Sample 16	Sample 20	Sample 24	Sample 28	Sample 32	Sample 36	Sample 40
<b>H</b>	<b>C3</b>	<b>C7</b>	Sample4	Sample8	Sample 12	Sample 16	Sample 20	Sample 24	Sample 28	Sample 32	Sample 36	Sample 40



LEGENDplex™ Kits are manufactured by **BioLegend**  
8999 BioLegend Way  
San Diego, CA 92121  
Tel: 1.858.768.5800  
Tel: US & Canada Toll-Free: 1.877.Bio-Legend (1.877.246.5343)  
Fax: 1.877.455.9587  
Email: [info@biolegend.com](mailto:info@biolegend.com)  
[biolegend.com](http://biolegend.com)

For a complete list of world-wide BioLegend offices and distributors,  
please visit our website at: [biolegend.com](http://biolegend.com)