


**Alexa Fluor® 594 anti-mouse Galectin-9 Antibody**

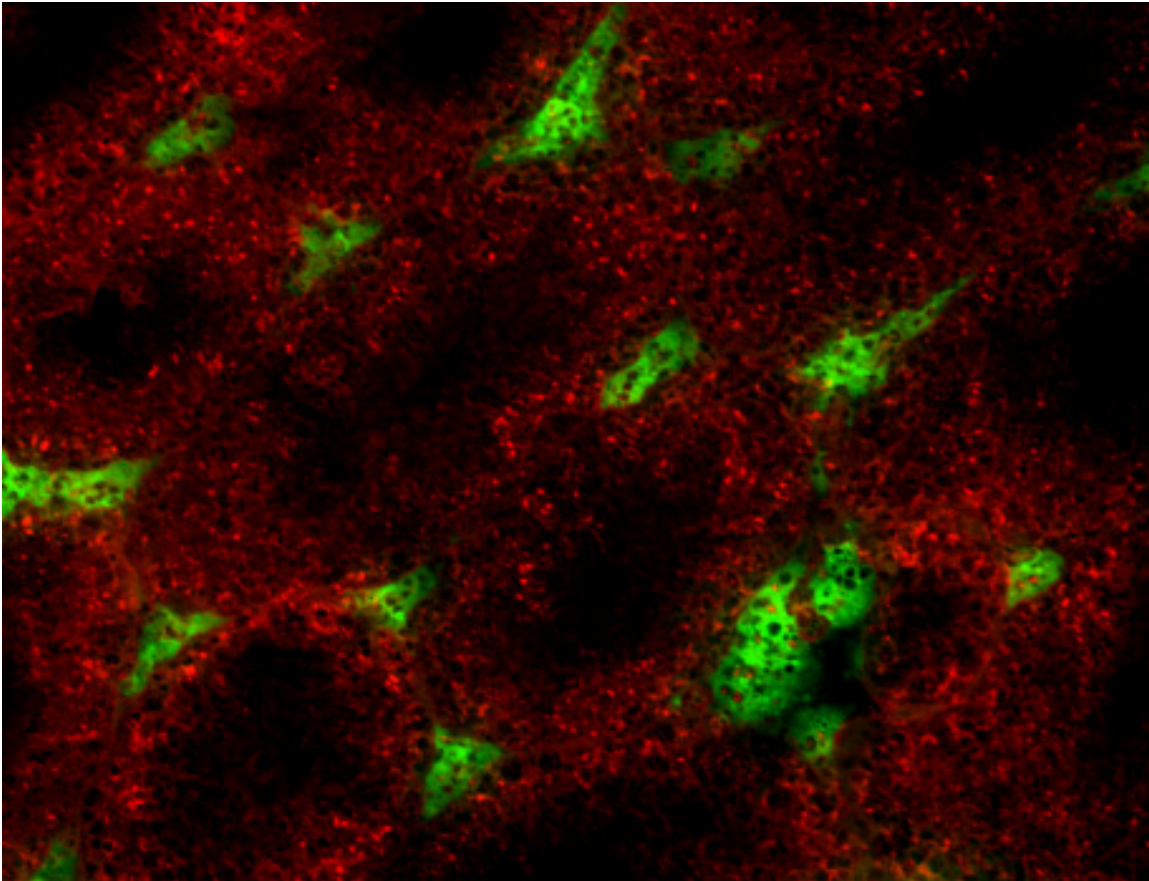

C57BL/6 mouse frozen testis section was fixed with 4% paraformaldehyde (PFA) for ten minutes at room temperature and blocked with 5% FBS plus 5% rat serum for one hour at room temperature. Then, the section was stained with 10 µg/ml of Galectin-9 (clone 108A2) Alexa Fluor® 594 (red) and 5 µg/mL of CD106 (clone 429 (MVCAM.A)) Alexa Fluor® 647(green) overnight at 4°C. The image was captured with a 10X objective.