


**Purified anti-Neurofilament H & M (NF-H/NF-M), Phosphorylated Antibody**

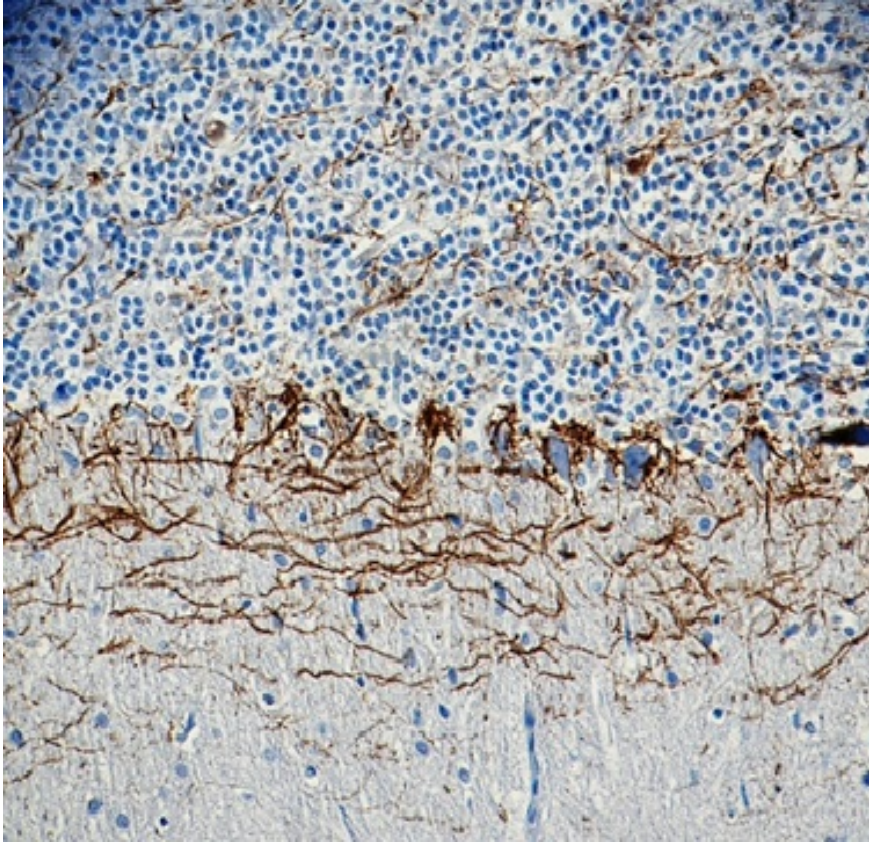

IHC staining of clone SMI 310 on formalin-fixed paraffin-embedded human cerebellum tissue at 0.5  $\mu\text{g}/\text{ml}$ . Following Retrieve-All Antigen Unmasking System 3: Acidic, 10X, the primary antibody was incubated for 60 minutes at room temperature. BioLegend's Ultra Streptavidin (USA) HRP Detection Kit (Multi-Species, DAB) was used for detection with DAB chromogen.